

# PHARMACIA

Pharmacia Groningen BV  
Van Swietenlaan 5  
9728 NX GRONINGEN  
The Netherlands

Phone Nr.: + 31 50 5296600  
Fax Nr.: + 31 50 5267860

## REPORT

TITLE: Bifocal foldable lens design based on  
corneal wavefront aberration

DATE: [REDACTED]

REPORT#: Rp2691r

AUTHOR: FUNCTION:  
Henk Weeber Research Scientist Biophysics

SIGNATURE: DATE:

[REDACTED]

APPROVED: FUNCTION:  
Marrie v.d. Mooren Sr. Development Engineer

SIGNATURE: DATE:

[REDACTED]

APPROVED: FUNCTION:  
Thom Terwee Manager Biophysics Department

SIGNATURE: DATE:

[REDACTED]

DISTRIBUTION: DOCUMENTATION CONTROL (CENTRAL FILE)  
Project file 513  
John Kappelhof  
Sverker Norrby

## TABLE OF CONTENTS

	Page
1. EXECUTIVE SUMMARY	3
2. INTRODUCTION	4
3. METHODS	4
3.1. Design methods	5
3.1.1. Design cornea	6
3.1.2. Light distribution	6
3.1.3. Power add for near vision	6
3.1.4. Base power	6
3.1.5. Optical quality	7
3.1.6. Optimization of optical quality	7
3.2. Measurement methods	8
3.2.1. Z9000 eye model	8
3.2.2. Light distribution	8
3.2.3. Power	8
3.2.4. Optical quality	9
4. RESULTS	9
4.1. Prototype designs	9
4.2. Theoretical performance of the prototype designs	10
4.3. Measurement of prototype lenses	13
4.3.1. Eye model	13
4.3.2. Modulation Transfer Function	15
4.3.3. Light distribution	16
4.3.4. Lens power	17
4.3.5. Lens power fine-tuning	17
5. DISCUSSION	18
6. REFERENCES:	19
Appendix 1 Design cornea	20
Appendix 2 OSLO Zernike operands	22
Appendix 3 PMMA Z9000 cornea in ISO fixture	24
Appendix 4 Measurement procedure of the lens power in the waterbath	30
Appendix 5 Test of the light distribution with different step heights	33
Appendix 6 Refractive indices	34
Appendix 7 Raw data	35
Appendix 8 CD rom: data files	56

## 1. EXECUTIVE SUMMARY

With Tecnis model Z9000, a lens is designed which corrects the corneal spherical aberration. The lens has an aspherical anterior surface, which induces a negative spherical aberration, equal but opposite to the average corneal spherical aberration of a cataract population. The resulting design has improved optical performance, especially for larger apertures<sup>1,2</sup>.

The design principle of Tecnis model Z9000 is applied on a bifocal silicone HRI lens. The principle bifocal specifications were copied from CeeOn model 811E: The design has a +4D reading add and a 50:50% light distribution between near and far vision. These specifications are met using a diffractive profile on the posterior lens surface.

For the posterior surface, the diameters of the rings are identical to CeeOn model 811E. The stepheight of the rings is 2.32 micrometer (model 811E: 1.85 micrometer).

Similar to Tecnis model Z9000, the anterior surface is used for the correction of spherical aberration. Two alternatives were evaluated:

1. an optimized anterior surface, which minimized the spherical aberration for near and far vision
2. an anterior surface identical to Tecnis model Z9000

The theoretical performance as well as measurements on prototype lenses revealed that there is no significant difference in performance between these two alternatives.

### Conclusions:

1. A bifocal foldable diffractive design with a regular Tecnis Z9000 anterior surface corrects for spherical aberration.
2. Compared to CeeOn model 811E, the resulting design shows approximately a factor 3 improvement of optical quality (MTF) for a 5mm pupil. For a 3mm pupil, the improvement is insignificant.

## 2. INTRODUCTION

With model Z9000, a lens is designed which corrects the corneal spherical aberration. The lens has an aspherical anterior surface, which induces a negative spherical aberration, equal but opposite to the average corneal spherical aberration of a cataract population. The resulting design has improved optical performance, especially for larger apertures<sup>1,2</sup>.

The design principle of model Z9000 can also be applied for bifocal lenses. Pharmacia has a PMMA bifocal lens, model 811E, which is a diffractive lens with a +4D power add for reading. This design can be transferred to HRI material<sup>3</sup>. Bifocal lenses have 2 foci. This means that one object will generate two images. In the ideal situation, one of the two images will be focused on the retina. While the other image is out of focus, it will still cause some blurr on the retina, which decreases the optical quality. Therefore, by their nature, the optical quality of bifocal lenses is lower than monofocal lenses of the same design.

This report shows how the design methods of model Z9000 are applied on a bifocal foldable lens. All files that were used to compile this report, including measurement data, calculation results and OSLO CCL programs, are stored on a CD-ROM as appendix 8 of this report.

## 3. METHODS

The design is based on the diffractive bifocal design of model 811E. This means that the same equations for the surface profile are used for the posterior surface, however with different parameters<sup>3</sup>(figure 1). Similar to model Z9000, the anterior surface is made aspherical<sup>4</sup>. Using these principles, the optical design of bifocal lenses comprises the determination of 3 items:

1. secure the light distribution between the 2 foci
2. secure the correct power add for near vision
3. secure the correct base power

With the new design, an additional item is added:

4. secure a specific spherical aberration

As for model 811E, the target light distribution between near and far focus is 50%:50% and the target power add for near vision is 4 diopters.

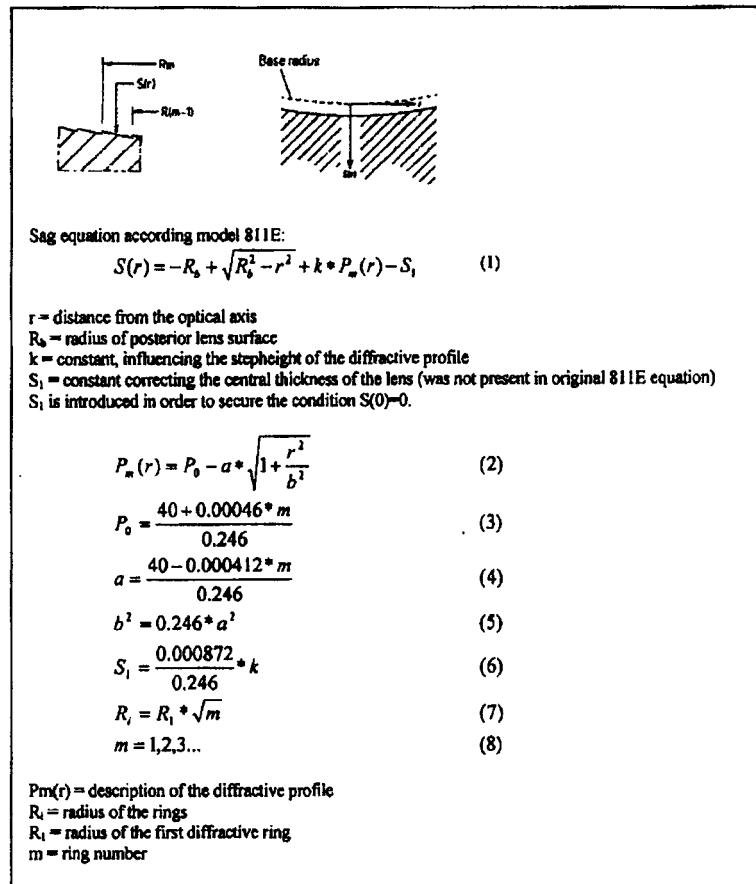


Figure 1. Equations describing the diffractive profile.

### 3.1. Design methods

The calculation methods which are used for determining the theoretical performance are described in detail in the report Rp2681r<sup>5</sup>.

### 3.1.1. Design cornea

The design cornea is described in report 2278. This design cornea is a 1-surface model. The refractive index of the cornea is the keratometry index of 1.3375. However, for diffractive lenses it is essential to use the real in vivo refractive index on the posterior (diffractive) lens surface. Therefore, a 2-surface model is used which has the same characteristics of the 1-surface model.

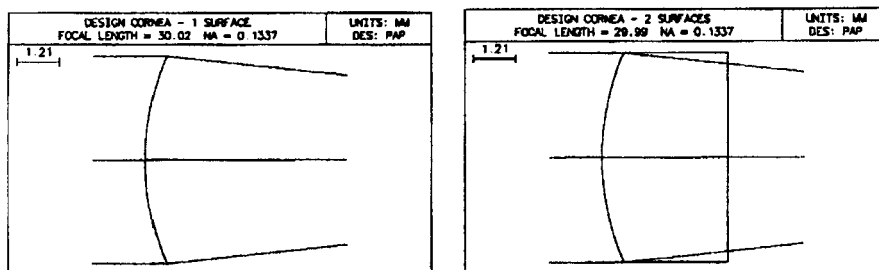


Figure 2. 1-surface (left) and 2-surface (right) design cornea.

The 2-surface cornea was constructed from the 1-surface cornea by adding a flat transition (surface) between the cornea and an aqueous/vitreous at 3.6mm from the cornea apex. The spherical aberration of the eye is remained within 0.03% of the original 1-surface cornea (see appendix 1).

### 3.1.2. Light distribution

The light distribution of the diffractive bifocal lens is determined by the stepheight of the diffractive zones.

The theoretical stepheight for a 50:50% light distribution is  $2.25 \mu\text{m}^3$ . This theoretical value is challenged and fine-tuned by means of prototype testing. It has to be noted that the light distribution is now defined in the Z9000 eye model (design cornea), while for model 811E, the light distribution was defined in the Gullstrand eye model.

### 3.1.3. Power add for near vision

The power add is determined by the diameters of the diffractive zones. Theoretically, this is independent of the refractive indices of the lens and the surrounding medium. Therefore, the same zone diameters are used as for model 811E.

### 3.1.4. Base power

The lens power of a normal monofocal lens is defined by the paraxial focus in situ<sup>6,7</sup>. Diffractive bifocal lenses do not have a theoretical paraxial focus. Therefore, an alternative approach was used for model 811E<sup>8</sup>. This approach basically means that, comparing a bifocal

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

---

and monofocal lens of the same lens power, the best focus should be located at the same position when the lens is measured in situ. This makes the practical behavior of monofocal and bifocal lens powers comparable.

In the design phase of the lens, an estimate can be generated by giving the diffractive lens a base radius comparable to the radius of a monofocal design. Fine-tuning is done in a later stage, when prototype lenses are measured.

### 3.1.5. Optical quality

The optical quality of diffractive bifocal lenses can be estimated using the optical design program OSLO<sup>9</sup>.

The methods are described in detail in report Rp2681r<sup>5</sup>.

### 3.1.6. Optimization of optical quality

The anterior surface of the bifocal lens can be optimized in order to reduce the overall spherical aberration, equivalent to what was done for model Z9000<sup>1</sup>. Alternatively, the exact Z9000 anterior surface can be applied onto the bifocal lens.

Optimization of the anterior surface is performed on the symmetrical Zernike terms:

OSLO Zernike term	weight factor
Z8	1
Z15	0.1
Z24	0.01

The optimization is performed for far and near vision simultaneously. For the calculations in OSLO, it means that for each optimization step, the vector sum of the Zernike term for near vision and far vision is taken as the operand (appendix 2). A real near vision configuration was generated by putting the light source in a position close to the eye.

As a starting point, the model Z9000 conic constant is used and  $ad=ae=0$ . The polynomial terms  $ad$  and  $ae$  are used for optimization.

### 3.2. Measurement methods

#### 3.2.1. Z9000 eye model

For the measurement of the optical modulation transfer function (MTF), a new cornea was used, which fits in front of a wet cell. The wet cell is identical to the ISO wet cell. The new PMMA cornea generates a wavefront that is identical to the wavefront of the design cornea. The focal point is in air, which facilitates automated focussing on the Eros Ealing optical bench. The design of the cornea is described in appendix 3. The refractive index of PMMA is measured from the actual PMMA rod. The refractive index of water is adjusted to match with the silicone HRI material: For diffractive lenses, the difference between the refractive index of the lens and the medium (aqueous or water) is important for the light distribution of the lens. Therefore, water with a refractive index of 1.3406 has to be used (Appendix 6). The refractive index of the water is changed by adding poly-ethyleenglycol (PEG200).

The advantages of the eye model over the previous one (report 2236) are:

- the PMMA cornea is not in direct contact to the medium (water), which makes it easier to maintain.
- The focus position is in air, so that the automated focusing of the Ealing optical bench can be used.

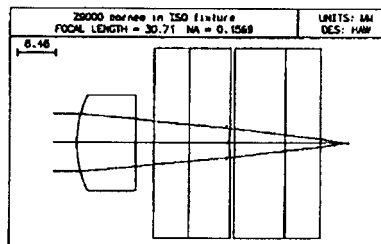


Figure 3. Z9000 cornea in front of the ISO wet cell.

#### 3.2.2. Light distribution

The light distribution is measured in the Z9000 eye model. At each focal point, the line-spread-function is measured. The central four pixel intensities are used as a measure of the amount of light in each focus.

#### 3.2.3. Power

The power of the bifocal lens is determined by measuring the back focal length (BFL) in the waterbath. The conversion of BFL to lens power is based on the following definition of bifocal lens power<sup>8</sup>:

The bifocal power of the lens is calculated from the back focal length in best focus condition, in the same way as is done for the paraxial lens power of a monofocal lens with equal geometry (radii and central thickness).

Calculation to be performed under the conditions:

1. The spherical aberration is presumed not to be influenced by diffractive pattern.
2. The location of the second principle plane is presumed not to be influenced by the diffractive pattern and is the same for both focal points.

Measurement of best focus to be determined under the conditions:



1. Collimated light (wavelength 550 nm) reaches the lens.
2. There is a 3mm aperture in front of the lens.
3. Best focus is defined as the position with the highest MTF-value at a spatial frequency of 50c/mm.

There is one difference with the original definition which was used for model 811E: For the new model, the spherical aberration depends on the lens power, and therefore on the focal point that is measured. In practice it means that for the far- and near focus different values for the spherical aberration have to be used. The background is that the spherical aberration of the current lens is much higher (when measured in a waterbath) than a spherical lens, while the amount of spherical aberration is also highly power-dependent. For model 811E, this diopter dependency was negligible.

The add power is defined as the difference between the powers of the near focus and far focus.

The measurement procedure is described in appendix 4.

### 3.2.4. Optical quality

The MTF is measured in the Z9000 eye model. Focussing is performed at 25c/mm, which is according to the current draft standards for multifocal lenses. The lenses are measured with apertures of 3 and 5 mm diameter. At 5mm aperture, the lens is measured at 3 positions, 45 degree apart, in order to account for any astigmatism.

## 4. RESULTS

### 4.1. Prototype designs

The prototype designs consist of an anterior and posterior optical surface. The posterior surface is the diffractive surface, according the equation also used for model 811E. The parameters  $k$  and  $R_1$  (equations 1 and 7 in figure 1), defining stepheight and the ring-diameters, have to be adapted for the silicone HRI design. Theoretically, the ring diameters for a given power add are independent of the material used, so the value used for model 811E is taken ( $R_1=0.513\text{mm}$ ).

The  $k$ -value determines the light distribution. In order to verify the  $k$ -value, a separate test series was made and tested (appendix 5). This test indicated that a  $k$ -value of 0.6536 results in a 50:50 light distribution. The corresponding stepheight of the rings is 2.32 micrometer.

The anterior surface was optimized for the symmetrical Zernike terms in far and near vision. The optimization has been performed for lens powers 15.0, 20.0 and 26.0D. On average, the resulting spherical aberration (OSLO Zernike Z8) is reduced with a factor 100. The resulting Zernike coefficients are listed in figure 4, the lens geometry in figure 5.

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

Power	Anterior surface	FAR			NEAR		
		Z8	Z15	Z24	Z8	Z15	Z24
15D	Optimized	-0.04	-0.11	-0.02	0.02	-0.10	-0.02
	Z9000	-0.07	0.01	-0.02	-0.01	0.01	-0.02
	Spherical	1.41	0.03	-0.02	1.46	0.03	-0.02
20D	Optimized	-0.03	-0.12	-0.02	0.02	-0.11	-0.02
	Z9000	-0.08	0.00	-0.02	-0.04	0.01	-0.02
	Spherical	1.66	0.04	-0.02	1.70	0.04	-0.02
26D	Optimized	-0.02	-0.12	-0.03	0.01	-0.12	-0.02
	Z9000	-0.09	0.00	-0.02	-0.06	0.01	-0.02
	Spherical	2.12	0.06	-0.02	2.16	0.06	-0.02

Figure 4. Theoretical performance of the prototype design, in terms of the symmetric Zernike coefficients. The Zernike coefficients are expressed in wavelengths. They are calculated in OSLO, in the design eye model, with a rectangular spot diagram in image space with using over 200.000 rays. For near vision, the reading distance is 40cm. For far vision, the aperture at the cornea is 6mm in diameter, which corresponds with an aperture diameter of 5.1mm at the lens position. For near vision, the entrance beam radius is adapted: The pupil size at the IOL location is the same for far and near vision.

	15D	20D	26D	
Anterior surface				
Radius	16.225	12.154	9.333	mm
Conic constant	-0.98856	-1.01855	-1.05018	H
AD	-4.3209E-04	-4.8285E-04	-6.0581E-04	mm <sup>3</sup>
AE	-4.2491E-05	-4.8990E-05	-5.9954E-05	mm <sup>4</sup>
Posterior surface				
Geometric radius	(-16.225)	(-12.154)	(-9.333)	mm
Base radius R <sub>0</sub>	-17.082	-12.586	-9.547	mm
Ring radius R <sub>1</sub>	0.513	0.513	0.513	mm
k-value	0.6536	0.6536	0.6536	mm
stepheight	2.32	2.32	2.32	µm
Central thickness	1.03	1.13	1.24	mm

Figure 5. Lens geometry of the prototype lenses\*. The geometric radius is added to indicate the overall geometry of the lens as if the posterior radius is a normal spherical surface, in this case equi-biconvex.

#### 4.2. Theoretical performance of the prototype designs

In the previous paragraph, the theoretical performance of the prototype designs was described in terms of the symmetrical Zernike coefficients. Figure 4 shows that the optimised design results in the smallest spherical aberration (Z8). Also the design with the original Z9000 shows very small values for the spherical aberration, though a little higher than the optimised design.

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

The modulation transfer function can be evaluated of the prototype designs, compared with a spherical diffractive design and with a spherical monofocal lens design. (figure 6, 7).

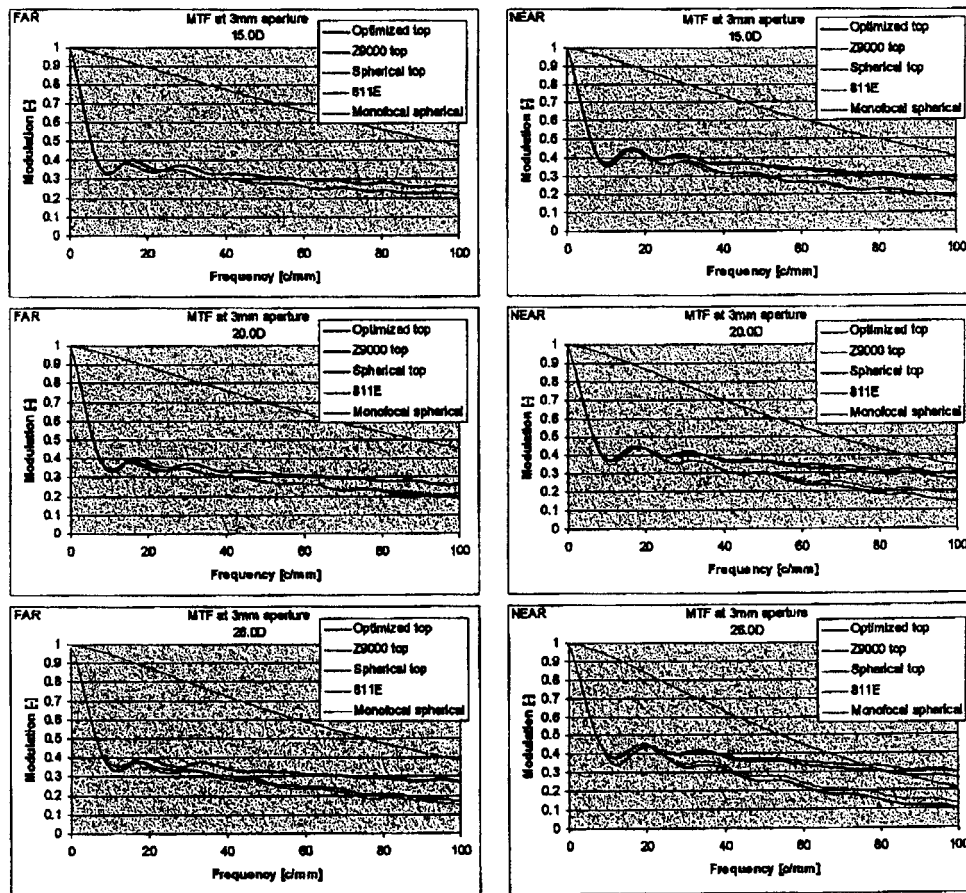


Figure 6. modulation transfer function of the prototype designs, compared with a spherical diffractive design and with a spherical monofocal lens design at 3mm.

The figures 6 and 7 show that the gain of the aspherical surface is there for the 5mm aperture, where the bifocal lens shows a contrast even higher than a spherical monofocal lens. Furthermore, the optimised anterior surface shows comparable MTF to the Z9000 anterior surface. For the very high frequencies at near vision the Z9000 surface is even slightly better.

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

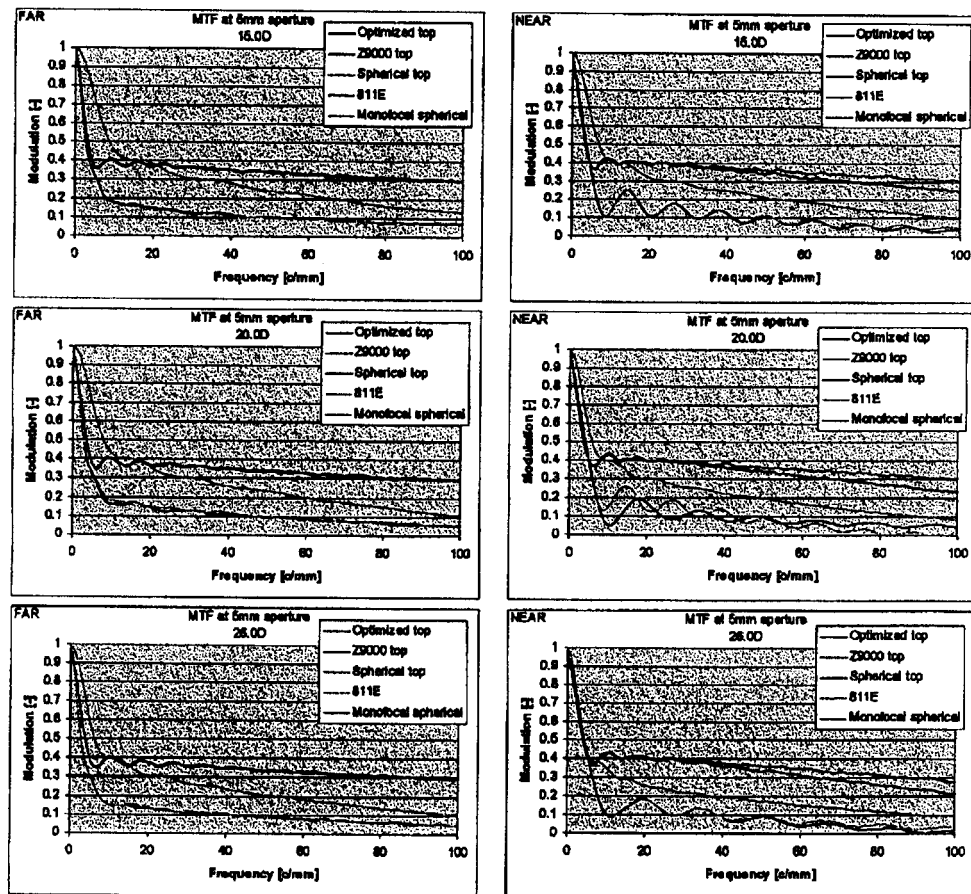


Figure 7. modulation transfer function of the prototype designs, compared with a spherical diffractive design and with a spherical monofocal lens design at 5mm. The heavy fluctuations of the MTF for the spherical top bifocal and 811E do not represent phase reversals of the phase transfer function. There is one phase reversal: model 811E, 20.0D, near vision, at 95c/mm (modulation = 0). Apart from this, the phase is always close to zero ( $<1^\circ$ ).

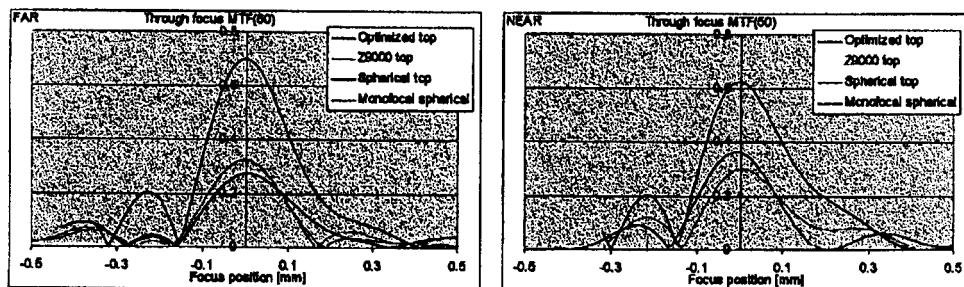


Figure 8. Through focus MTF at 50c/mm for 20.0D lenses and a 3mm aperture.

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

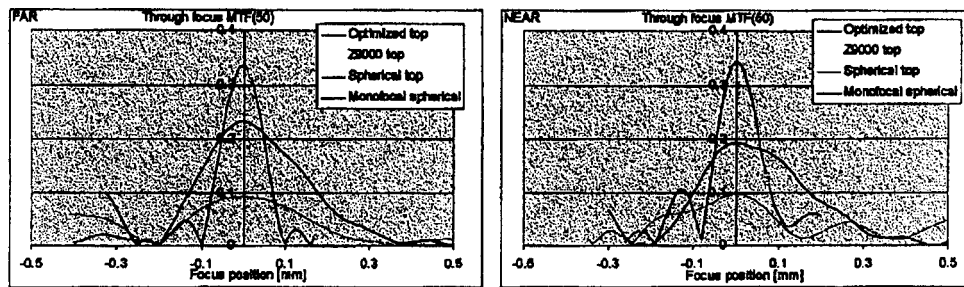


Figure 9. Through focus MTF at 50c/mm for 20.0D lenses and a 5mm aperture.

### 4.3. Measurement of prototype lenses

The lenses are measured on their optical performance, in an eye model and in a waterbath. Shape measurements on diffractive and aspherical surfaces are extremely difficult, especially on a flexible and sticky silicone surface. Feasibility measurements were carried out with optical and mechanical measuring techniques with varying success. The results were not reliable enough for extensive reporting here. Nevertheless all results showed that the lenses were close to the design in the sub-micron range (average deviation of stepheight 0.04 micrometer). Also mechanical tests on the molds (Talysurf) indicated a good result. Furthermore, the shape of the lenses was secured by using production procedures standard for regular Z9000 lenses. The test lenses were not angle set and not sterilized.

#### 4.3.1. Eye model

The lenses were measured in the new Z9000 eye model as described in appendix 3. The water/PEG200 solution with a refractive index of 1.3406 was made and measured in the AR chemical lab (appendix 6).

The Z9000 PMMA cornea in front of the ISO wet cell is qualified by through focus and through frequency MTF measurements at 3mm and 5mm pupil. The measurement results are displayed in figure 10, together with the theoretical curves. Focussing was performed at 50c/mm. The measured graphs compare well with the calculated graphs. As a check, a Z9000 lens was measured as well (figure 11).

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

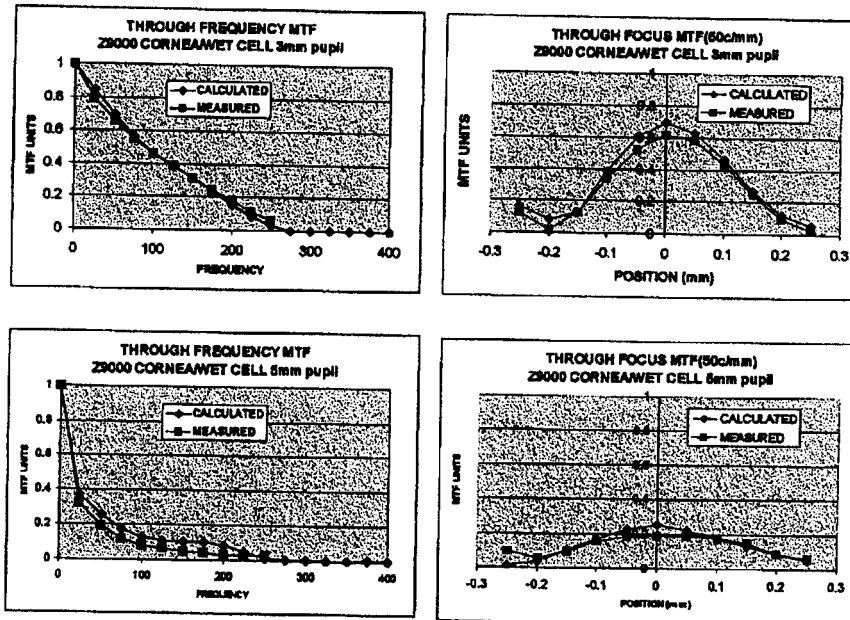


Figure 10. Measured through frequency and through focus MTF curves for the new Z9000 cornea for the wet cell, both for a 3mm and 5mm aperture.

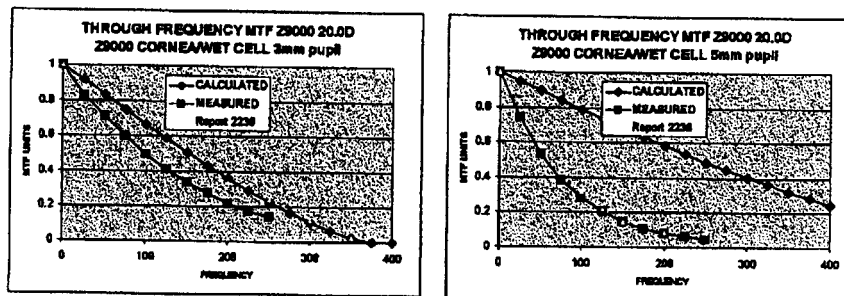


Figure 11. Measured through frequency MTF curves for a Z9000 lens, 20.0D, in the eye model. As a reference, the measured data from the lens verification report 2236 is added.

#### 4.3.2. Modulation Transfer Function

Lenses were measured in the eye model on the optical bench. The measuring conditions were:

- Refractive index of the water (+PEG200) 1.3406
- Automatic focussing on MTF at 25c/mm
- Maximum analysis window of the linear photodiode array, leaving  $1/16$  unused at the edges.
- For 3mm apertures, 1 measurement per focal point was made, for 5mm apertures, 3 measurements per focal point were made (lens at  $-45^\circ$ ,  $0^\circ$  and  $+45^\circ$ )
- 8 lenses per prototype, per lens power

The results are given in figures 12. The results show that the optimized top and the Z9000 top have comparable performance in terms of MTF. At a 5mm aperture they are much better than model 811E.

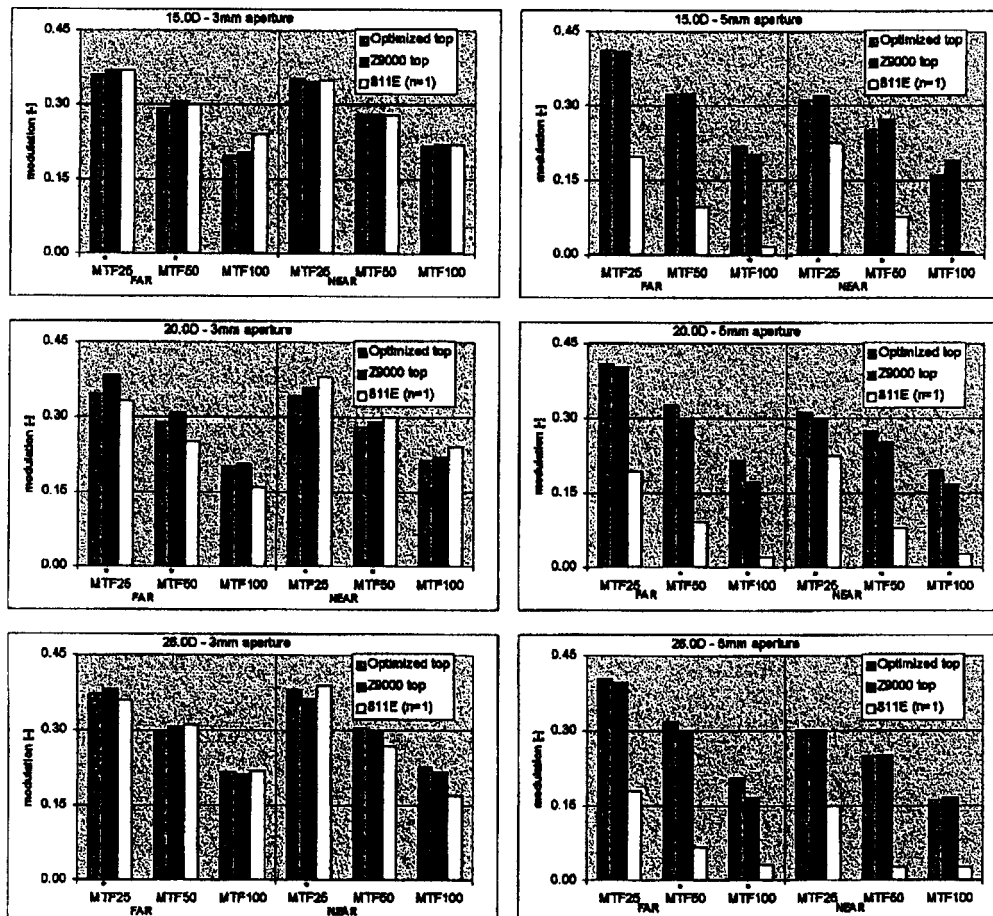


Figure 12. Modulation transfer values for 25, 50 and 100c/mm and aperture 3mm. The data of the optimized top and the Z9000 top is the average of 8 lenses. The data of model 811E is based on the measurement of one lens per diopter.

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

The results show that the optimized top and the Z9000 top have comparable performance in terms of MTF. At a 5mm aperture they are much better than model 811E. The standard deviation of every series of 8 lenses was small, with a standard deviation less than 0.02 MTF units (3mm and 5mm aperture). For the 5mm aperture, 3 measurements were performed per focal point, while the lens was rotated 45 degrees in between the measurements. The variation between these 3 measurements on the same lens depends on the spatial frequency and was 0.01, 0.02 and 0.03 MTF units (s.d.) for 25, 50 and 100c/mm respectively.

Figure 13 shows the ratio between the measured MTF and the theoretical MTF for the prototype with the Z9000 top. The theoretical MTF was calculated for lens under MO-lab conditions. The average ratio of the 3 spatial frequencies for the 3mm pupil is 88%, with a minimum of 71%. This is just within the ISO standard for monofocal lenses, which states that the lens optical quality should be greater or equal than 70% of the theoretical performance.

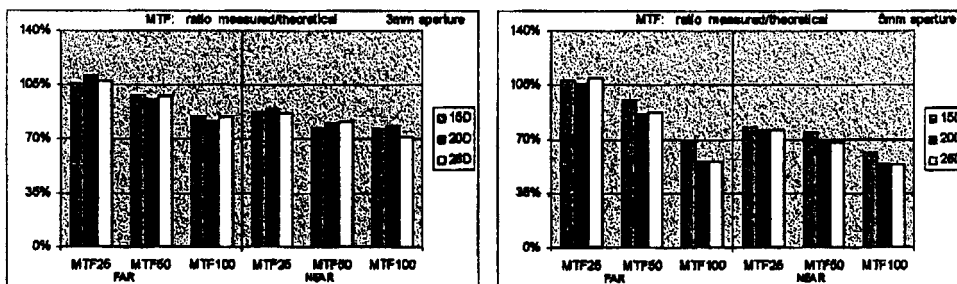


Figure 13. Ratio of measured/calculated MTF for prototype with Z9000 top.

#### 4.3.3. Light distribution

The light distribution was calculated from the signal on the diode array, obtained during the MTF measurements for a 3mm aperture. The light distribution is defined as the ratio between the far and near focus. Figure 14 shows only the amount going to the far focus. The data shows no relationship between light distribution and lens power.

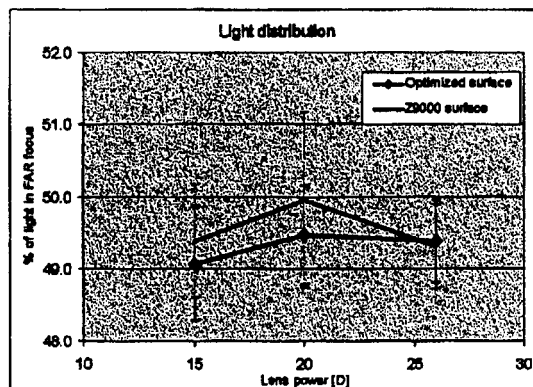


Figure 14. Light distribution: Amount of light that reaches the far focus. The error bars correspond with plus and minus one standard deviation.



#### 4.3.4. Lens power

The lens power was determined by measuring 8 lenses per diopter. Only the prototype with the Z9000 top has been measured. The add power is very well on the target of 4 diopter. The base power (far) is a little off and needs fine-tuning.

Nominal power	Measured lens power			Deviation	
	FAR	NEAR	ADD	FAR	ADD
15.0	14.55	18.53	3.98	-0.45	-0.02
20.0	19.66	23.71	4.05	-0.34	0.05
26.0	25.78	29.75	3.97	-0.22	-0.03

Nominal power	Standard deviations [D]		
	FAR	NEAR	ADD
15.0	0.08	0.14	0.15
20.0	0.11	0.04	0.10
26.0	0.05	0.06	0.03

Figure 15. Measured lens power and standard deviations.

#### 4.3.5. Lens power fine-tuning

Based on the measurement results, the lens power needs fine-tuning. The method used here, is to compare the measured result with a calculated paraxial lens power, based on the design anterior radius, the design central thickness and the design base posterior radius  $R_b$  (equation 1, figure 1). This leads to the following deviation of lens powers:

Nominal Power	design values			TLE* Power	Measured power	deviation	proposed correction
	Anterior radius	Posterior base radius $R_b$	Central Thickness				
15.0	18.225	-17.082	1.03	14.62	14.55	-0.07	0.07
20.0	12.154	-12.586	1.13	19.66	19.66	0.00	0.00
26.0	9.333	-9.547	1.24	25.71	25.78	0.07	-0.07

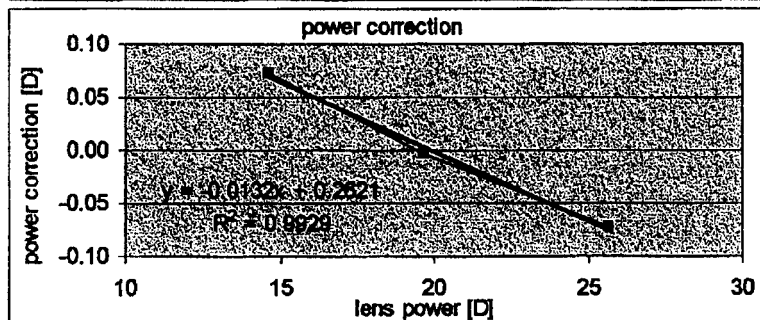


Figure 16. Lens power correction of the prototype with Z9000 top.

\* TLE = thick-lens-equation, using the design geometry as input data.

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

---

The lens power correction is applied on the base radius Rb. The resulting lens data are:

Lens Power	Anterior surface identical to Z9000				Posterior surface k=0.6536 R1=0.513	
	Radius	CC	AD	AE	Base radius Rb	Central Thickness
15.0	16.225	-0.988562	-8.03E-04	-7.23E-06	-16.087	1.03
15.5	15.899	-0.985816	-8.12E-04	-7.57E-06	-15.583	1.04
16.0	15.207	-0.99036	-8.21E-04	-7.93E-06	-15.110	1.05
16.5	14.745	-0.99753	-8.29E-04	-8.30E-06	-14.665	1.06
17.0	14.309	-0.99846	-8.39E-04	-8.69E-06	-14.245	1.07
17.5	13.899	-1.00088	-8.48E-04	-9.09E-06	-13.848	1.08
18.0	13.511	-1.00077	-8.58E-04	-9.52E-06	-13.473	1.09
18.5	13.144	-1.00030	-8.68E-04	-9.96E-06	-13.117	1.10
19.0	12.797	-1.00744	-8.78E-04	-1.04E-05	-12.780	1.11
19.5	12.467	-1.00314	-8.89E-04	-1.09E-05	-12.459	1.12
20.0	12.154	-1.018548	-8.99E-04	-1.14E-05	-12.154	1.13
20.5	11.856	-1.023293	-9.10E-04	-1.19E-05	-11.864	1.13
21.0	11.572	-1.027614	-9.21E-04	-1.25E-05	-11.587	1.14
21.5	11.301	-1.030542	-9.33E-04	-1.31E-05	-11.322	1.15
22.0	11.043	-1.036125	-9.44E-04	-1.37E-05	-11.069	1.16
22.5	10.796	-1.034876	-9.57E-04	-1.43E-05	-10.828	1.17
23.0	10.560	-1.037706	-9.69E-04	-1.50E-05	-10.596	1.18
23.5	10.334	-1.039948	-9.81E-04	-1.56E-05	-10.374	1.19
24.0	10.117	-1.041622	-9.94E-04	-1.63E-05	-10.161	1.20
24.5	9.909	-1.043677	-1.01E-03	-1.71E-05	-9.957	1.21
25.0	9.709	-1.046434	-1.02E-03	-1.79E-05	-9.760	1.22
25.5	9.517	-1.048667	-1.04E-03	-1.87E-05	-9.572	1.23
26.0	9.333	-1.050179	-1.05E-03	-1.95E-05	-9.390	1.24

Figure 17. Lens optical design of the bifocal foldable lens with an anterior surface profile which is identical to Tecnis model Z9000.

## 5. DISCUSSION

This report has shown that the Tecnis Z9000 design principles can be successfully applied on bifocal lenses. Two approaches were used: one using the proven Z9000 anterior lens shape combined with a diffractive posterior surface. Alternatively a new anterior lens shape was generated by optimizing the wavefront aberrations for the far as well as the near focus. The performance of these two types of lenses, in terms of MTF, showed to be identical, in theory as well as according the measurement of prototype lenses.

The improvement of the ZM001, compared to model 811E, is significant. However this is only true for the larger pupils (larger than 3mm). For a 5mm aperture, the MTF of the new bifocal design is even better than a spherical monofocal lens. Caution has to be taken to

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

---

translate this to clinical behavior, since the clinical behavior is also influenced by:

- perception of the out-of-focus image
- scatter

The theoretical through focus behavior (figures 8 and 9) indicates that the depth of focus is somewhat reduced for the larger pupil sizes. However, compared to the monofocal Z9000 model there doesn't seem to much difference<sup>2</sup>. Through focus measurements of prototype lenses can give more insight in the through focus behavior. It should also be noted that the theoretical calculations of MTF of diffractive multifocal lenses have to be interpreted with caution. The calculation routines of out-of-focus MTF have not been verified. For similar reasons, no calculations have been performed with tilts and decentrations.

Currently, there are no approved standards for the optical performance of multifocal standards. For model 811E, Pharmacia has generated its own specification, based on measurement in the Gullstrand eye model. A similar approach could be performed with the new model, based on the Z9000 eye model. For the 3mm aperture, the measured MTF is more than 70% of its predicted behavior. In this respect it fulfils the ISO standard for monofocal lenses. Whether this is an appropriate specification would need to be further verification, since there is currently not much room left (71% at 26D). Deviation from this standard can be justified by:

1. The standard is for monofocal lenses, not bifocal lenses
2. Prediction of diffractive bifocal lenses is difficult, with limited accuracy
3. The Z9000 eye model has an aspherical surface, creating an additional measurement error

## 6. REFERENCES:

1. P. Piers, 'Z11 design based on corneal wavefront aberration', PG report R2278.
2. M. van der Mooren, 'CeeOn 911Z verification report', PG report Rp2236r.
3. H. Weeber, 'Feasibility report bifocal silicone IOL', PG report 1727.
4. H. Weeber, 'Influence of the position of the aspheric surface on optical quality', PG report Rp2524r
5. H. Weeber, 'Design tools for diffractive bifocal lenses', PG report Rp2681r.
6. ISO standard: ISO11979-1
7. ANSI standard: ANSI Z80.7-1994
8. J. Hermans, 'Bifocal lens power: Definition and practical implementation', PG report 1031.
9. OSLO Premium Edition Revision 6.1. Lambda Research Corporation, Littleton, MA

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

Appendix 1: Design cornea

\*LENS DATA

DESIGN CORNEA - 1 SURFACE

SRF	RADIUS	THICKNESS	APERTURE RADIUS	GLASS	SPE	NOTE
OBJ	--	1.0000e+20	1.0000e+14	AIR		
1	7.575000	3.600000	3.000003 S	KERATOM	M *	
AST	--	--	2.640233 AS	KERATOM	P	
3	--	0.900000	2.640233 S	KERATOM	P	
4	--	25.519444 S	2.550292 S	KERATOM	P	
IMS	--	--	2.2444e-05 S			

\*CONIC AND POLYNOMIAL ASPHERIC DATA

SRF	CC	AD	AE	AF	AG
1	-0.141350	--	--	--	--

\*REFRACTIVE INDICES

SRF	GLASS	RN1	RN2	RN3	VNBR	TCE
0	AIR	1.000000	1.000000	1.000000	--	--
1	KERATOM	1.337500	1.676807	1.207932	0.719807	236.000000
2	KERATOM	1.337500	1.676807	1.207932	0.719807	236.000000
3	KERATOM	1.337500	1.676807	1.207932	0.719807	236.000000
4	KERATOM	1.337500	1.676807	1.207932	0.719807	236.000000
5	IMAGE SURFACE					

\*WAVELENGTHS

CURRENT	WV1/WW1	WV2/WW2	WV3/WW3
1	0.587560	0.486130	0.656270
	1.000000	1.000000	1.000000

\*LENS DATA

DESIGN CORNEA - 2 SURFACES

SRF	RADIUS	THICKNESS	APERTURE RADIUS	GLASS	SPE	NOTE
OBJ	--	1.0000e+20	1.0000e+14	AIR		
1	7.575000	3.600000	3.000003 S	KERATOM	M *	
AST	--	--	2.640233 AS	KERATOM	P	
3	--	0.900000	2.640233 S	AQ	M	
4	--	25.489815 S	2.550191 S	VIT	M	
IMS	--	--	2.2444e-05 S			

\*CONIC AND POLYNOMIAL ASPHERIC DATA

SRF	CC	AD	AE	AF	AG
1	-0.141350	--	--	--	--

\*REFRACTIVE INDICES

SRF	GLASS	RN1	RN2	RN3	VNBR	TCE
0	AIR	1.000000	1.000000	1.000000	--	--
1	KERATOM	1.337500	1.676807	1.207932	0.719807	236.000000
2	KERATOM	1.337500	1.676807	1.207932	0.719807	236.000000
3	AQ	1.336000	1.675279	1.206444	0.716671	236.000000
4	VIT	1.336000	1.675279	1.206444	0.716671	236.000000
5	IMAGE SURFACE					

\*WAVELENGTHS

CURRENT	WV1/WW1	WV2/WW2	WV3/WW3
1	0.587560	0.486130	0.656270
	1.000000	1.000000	1.000000

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

---

The Z11 cornea is a 1-surface eye model. For the lens design, a 2-surface model is used in order to have the lens in aqueous. It was verified that the 2 eye models have the same wavefront aberrations:

(OSLO notation)

coefficient	1 surface	2 surface	Difference
Z0	0.001173714	0.001174047	0.03%
Z1	0	0	
Z2	0	0	
Z3	0.001769161	0.001769673	0.03%
Z4	0	0	
Z5	0	0	
Z6	0	0	
Z7	0	0	
Z8	0.000801207	0.000801396	0.03%
Z9	0	0	
Z10	0	0	
Z11	0	0	
Z12	0	0	
Z13	0	0	
Z14	0	0	
Z15	5.78493E-08	5.79E-08	0.16%
Z16	0	0	
Z17	0	0	
Z18	0	0	
Z19	0	0	
Z20	0	0	
Z21	0	0	
Z22	0	0	
Z23	0	0	
Z24	2.40715E-08	2.43E-08	0.76%
Z25	0	0	
Z26	0	0	
Z27	0	0	
Z28	0	0	
Z29	0	0	
Z30	0	0	
Z31	0	0	
Z32	0	0	
Z33	0	0	
Z34	0	0	
Z35	-5.60344E-10	-5.48E-10	-0.41%
Z36	0	0	

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

---

## Appendix 2: OSLO Zernike operands

The operands, used for optimizing the lens are calculated in a procedure 'oprd\_zernike2'.  
The procedure is stored in the file 'optimize\_Z\_haw.ccl':

```
cmd oprd_zernike2()
// Define Zernike terms as operands for optimization in 2 configurations
(
  set_preference(output_text, off);
  cfg 1; // Configuration 1
  sop 0;tra all 1; // om te voorkomen dat OSLO vastloopt bij diffractielenzen
  spd mon none +SDAD 1; // SDAD=APDIV in lensfile (setup) bepalend voor de
  rekensnelheid!!
  ssbuf_reset();zer bf 1 36;
  Ocm[0] = a1**2; // OSLO Zernike Piston term Z0
  Ocm[1] = a2**2; // Z1 TILT
  Ocm[2] = a3**2; // Z2 TILT
  Ocm[3] = a4**2; // Z3 DEFOCUS
  Ocm[4] = a5**2; // Z4 ASTIGMATISM
  Ocm[5] = a6**2; // Z5 ASTIGMATISM
  Ocm[6] = a7**2; // Z6 COMA
  Ocm[7] = a8**2; // Z7 COMA
  Ocm[8] = a9**2; // Z8 3rd ORDER SPHERICAL ABERRATION
  Ocm[9] = a11**2; // Z9 TREFOIL
  Ocm[10] = a11**2; // Z10 TREFOIL
  Ocm[11] = a12**2; // Z11
  Ocm[12] = a13**2; // Z12
  Ocm[13] = a14**2; // Z13
  Ocm[14] = a15**2; // Z14
  Ocm[15] = a16**2; // Z15
  Ocm[16] = a17**2; // Z16
  Ocm[17] = a18**2; // Z17
  Ocm[18] = a19**2; // Z18
  Ocm[19] = a20**2; // Z19
  Ocm[20] = a21**2; // Z20
  Ocm[21] = a22**2; // Z21
  Ocm[22] = a23**2; // Z22
  Ocm[23] = a24**2; // Z23
  Ocm[24] = a25**2; // Z24
  Ocm[25] = a26**2; // Z25
  Ocm[26] = a27**2; // Z26
  Ocm[27] = a28**2; // Z27
  Ocm[28] = a29**2; // Z28
  Ocm[29] = a30**2; // Z29
  Ocm[30] = a31**2; // Z30
  Ocm[31] = a32**2; // Z31
  Ocm[32] = a33**2; // Z32
  Ocm[33] = a34**2; // Z33
  Ocm[34] = a35**2; // Z34
  Ocm[35] = a36**2; // Z35
  Ocm[36] = a37**2; // Z36
  set_preference(output_text, on);
  set_preference(output_text, off);
  cfg 2; // Configuration 2
  sop 0;tra all 1; // om te voorkomen dat OSLO vastloopt bij diffractielenzen
  spd mon none +SDAD 1; // SDAD=APDIV in lensfile (setup) bepalend voor de
  rekensnelheid!!
  ssbuf_reset();zer bf 1 36;
  Ocm[0] += a1**2; // OSLO Zernike Piston term Z0
  Ocm[1] += a2**2; // Z1 TILT
  Ocm[2] += a3**2; // Z2 TILT
  Ocm[3] += a4**2; // Z3 DEFOCUS
  Ocm[4] += a5**2; // Z4 ASTIGMATISM
  Ocm[5] += a6**2; // Z5 ASTIGMATISM
  Ocm[6] += a7**2; // Z6 COMA
  Ocm[7] += a8**2; // Z7 COMA
  Ocm[8] += a9**2; // Z8 3rd ORDER SPHERICAL ABERRATION
  Ocm[9] += a11**2; // Z9 TREFOIL
  Ocm[10] += a11**2; // Z10 TREFOIL
  Ocm[11] += a12**2; // Z11
  Ocm[12] += a13**2; // Z12
  Ocm[13] += a14**2; // Z13
```

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

---

```
Ocm[14] += a15**2; // Z14
Ocm[15] += a16**2; // Z15
Ocm[16] += a17**2; // Z16
Ocm[17] += a18**2; // Z17
Ocm[18] += a19**2; // Z18
Ocm[19] += a20**2; // Z19
Ocm[20] += a21**2; // Z20
Ocm[21] += a22**2; // Z21
Ocm[22] += a23**2; // Z22
Ocm[23] += a24**2; // Z23
Ocm[24] += a25**2; // Z24
Ocm[25] += a26**2; // Z25
Ocm[26] += a27**2; // Z26
Ocm[27] += a28**2; // Z27
Ocm[28] += a29**2; // Z28
Ocm[29] += a30**2; // Z29
Ocm[30] += a31**2; // Z30
Ocm[31] += a32**2; // Z31
Ocm[32] += a33**2; // Z32
Ocm[33] += a34**2; // Z33
Ocm[34] += a35**2; // Z34
Ocm[35] += a36**2; // Z35
Ocm[36] += a37**2; // Z36
set_preference(output_text, on);
}
```

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

**Appendix 3: PMMA Z9000 cornea in ISO fixture**

ISO uses a cornea in front of a wet cell. The ISO cornea can be replaced by a PMMA cornea. The PMMA cornea has to generate a wavefront equivalent to the wavefront in the Z9000 model eye. This is characterised by the working f-number and the wavefront aberrations (Zernike terms).

The Z9000 model eye is characterized by:

Aperture diameter of the cornea = 6.0 mm


Corresponding aperture radius at the IOL position: 2.640233 mm

Working f-number = 3.740741

Wavefront aberration (OSLO notation):  
Z8 = 0.000601 mm  
Z15 =  $5.7840 \times 10^{-6}$  mm  
(Z24 =  $2.4 \times 10^{-8}$  mm)


The ISO eyemodel is described in ISO/DIS 11979-2. Only the wet cell with the aperture is used, with pure water as medium. The PMMA cornea has a central thickness of 10.0mm and is placed 3.0mm in front of the cell. The refractive index of the rod PMMA was measured as 1.4872 (22°C, 546.1nm). For the calculation of the working f-number and the wavefront aberrations, the right window is taken away, so that the system focusses in water. The working f-number determines the radius of the convex-plano PMMA cornea:  $R=18.5683$  mm. The aspherical components of the PMMA surface are optimized to produce the Z9000 wavefront aberrations, which resulted in:

cc = -0.177999  
ad =  $+5.1185 \times 10^{-6}$   
ac =  $+3.3540 \times 10^{-9}$

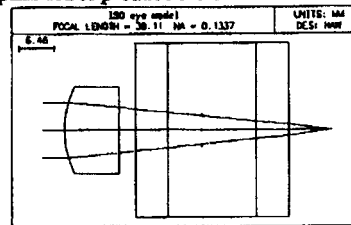

When the measurement is taken with the focus in air (the real situation), the resulting wavefront aberration is slightly higher (4% for Z8). When a 20.0 diopter Z9000 IOL is placed in the cell, 94% of the wavefront aberration (Z8) is corrected. In the Z9000 design eyemodel, a 20.0D IOL corrects 99% of the wavefront aberration (Z8).

MTF for a 5.1mm aperture (at the IOL position):

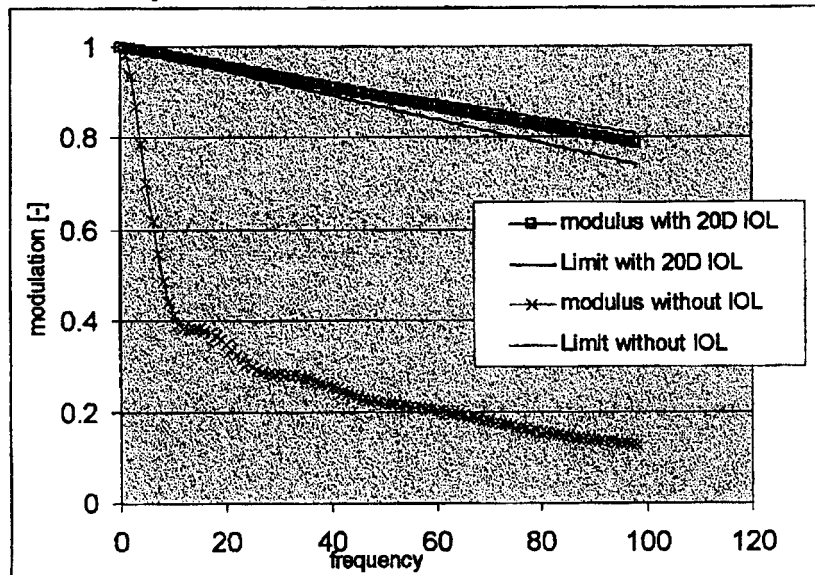



REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

---

**Verification: PMMA cornea for ISO-fixture, with focus in water.**

Apdiv 97.10  
Zernike in mm  
Spot diagram in image space  
Pupil 2.640233

**\*LENS DATA**

PMMA Z9000 cornea in ISO fixture

SRF	RADIUS	THICKNESS	APERTURE RADIUS	GLASS	SPE	NOTE
OBJ	--	1.0000e+20	1.0000e+14	AIR		
1	18.568300	10.000000	8.000000	PMMA_ROD	*	
2	--	3.000000	8.000000	AIR		
3	--	6.000000	16.000000	BK7	C	
4	--	6.250000	16.000000	WATER	C	
AST	--	26.355802 S	2.640233 AS	WATER	C	
IMS	--	--	3.8112e-05 S			

**\*CONIC AND POLYNOMIAL ASPHERIC DATA**

SRF	CC	AD	AE	AF	AG
1	-0.177999	5.1185e-06	3.3540e-09	--	--

**\*ZERNIKE ANALYSIS**

WAVELENGTH 1

Positive angle A is a rotation from the +y axis toward the +x axis.

0.001173: [0] 1  
-- : [1] RCOSA  
-- : [2] RSINA  
0.001769: [3] 2R^2 - 1  
-- : [4] R^2COS2A  
-- : [5] R^2SIN2A  
-- : [6] (3R^2 - 2)RCOSA  
-- : [7] (3R^2 - 2)RSINA  
0.000601: [8] 6R^4 - 6R^2 + 1  
-- : [9] R^3COS3A  
-- : [10] R^3SIN3A  
-- : [11] (4R^2 - 3)R^2COS2A  
-- : [12] (4R^2 - 3)R^2SIN2A  
-- : [13] (10R^4 - 12R^2 + 3)RCOSA  
-- : [14] (10R^4 - 12R^2 + 3)RSINA  
5.7840e-06: [15] 20R^6 - 30R^4 + 12R^2 - 1  
-- : [16] R^4COS4A  
-- : [17] R^4SIN4A  
-- : [18] (5R^2 - 4)R^3COS3A  
-- : [19] (5R^2 - 4)R^3SIN3A  
-- : [20] (15R^4 - 20R^2 + 6)R^2COS2A  
-- : [21] (15R^4 - 20R^2 + 6)R^2SIN2A  
-- : [22] (35R^6 - 60R^4 + 30R^2 - 4)RCOSA  
-- : [23] (35R^6 - 60R^4 + 30R^2 - 4)RSINA  
2.3253e-08: [24] 70R^8 - 140R^6 + 90R^4 - 20R^2 + 1  
-- : [25] R^5COS5A  
-- : [26] R^5SIN5A  
-- : [27] (6R^2 - 5)R^4COS4A  
-- : [28] (6R^2 - 5)R^4SIN4A  
-- : [29] (21R^4 - 30R^2 + 10)R^3COS3A  
-- : [30] (21R^4 - 30R^2 + 10)R^3SIN3A  
-- : [31] (56R^6 - 105R^4 + 60R^2 - 10)R^2COS2A  
-- : [32] (56R^6 - 105R^4 + 60R^2 - 10)R^2SIN2A  
-- : [33] (126R^8 - 280R^6 + 210R^4 - 60R^2 + 5)RCOSA  
-- : [34] (126R^8 - 280R^6 + 210R^4 - 60R^2 + 5)RSINA  
-6.0838e-10: [35] 252R^10 - 630R^8 + 560R^6 - 210R^4 + 30R^2 - 1  
-- : [36] 924R^12 - 2772R^10 + 3150R^8 - 1680R^6 + 420R^4 - 42R^2 + 1  
RMS OPD 0.001056 ERROR 7.0480e-11

**Verification: PMMA cornea for ISO-fixture, with focus in air.**

**\*LENS DATA**

PMMA Z9000 cornea in ISO fixture

SRF	RADIUS	THICKNESS	APERTURE RADIUS	GLASS	SPE	NOTE
-----	--------	-----------	-----------------	-------	-----	------

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

OBJ	--	1.0000e+20	1.0000e+14	AIR
1	18.568300	10.000000	8.000000	PMMA_ROD *
2	--	3.000000	8.000000	AIR
3	--	6.000000	16.000000	BK7 C
4	--	6.250000	16.000000	WATER C
AST	--	--	2.640233 AS	WATER C
6	--	10.000000	2.640233 S	WATER C
7	--	6.000000	16.000000	BK7 C
8	--	8.306948 S	16.000000	AIR
IMS	--	--	3.8112e-05 S	

\*CONIC AND POLYNOMIAL ASPHERIC DATA

SRF	CC	AD	AE	AF	AG
1	-0.177999	5.1185e-06	3.3540e-09	--	--

\*ZERNIKE ANALYSIS

WAVELENGTH 1

Positive angle A is a rotation from the +y axis toward the +x axis.

```

0.001212: [0] 1
-- : [1] RCOSA
-- : [2] RSINA
0.001828: [3] 2R^2 - 1
-- : [4] R^2COS2A
-- : [5] R^2SIN2A
-- : [6] (3R^2 - 2)RCOSA
-- : [7] (3R^2 - 2)RSINA
0.000623: [8] 6R^4 - 6R^2 + 1
-- : [9] R^3COS3A
-- : [10] R^3SIN3A
-- : [11] (4R^2 - 3)R^2COS2A
-- : [12] (4R^2 - 3)R^2SIN2A
-- : [13] (10R^4 - 12R^2 + 3)RCOSA
-- : [14] (10R^4 - 12R^2 + 3)RSINA
6.7871e-06: [15] 20R^6 - 30R^4 + 12R^2 - 1
-- : [16] R^4COS4A
-- : [17] R^4SIN4A
-- : [18] (5R^2 - 4)R^3COS3A
-- : [19] (5R^2 - 4)R^3SIN3A
-- : [20] (15R^4 - 20R^2 + 6)R^2COS2A
-- : [21] (15R^4 - 20R^2 + 6)R^2SIN2A
-- : [22] (35R^6 - 60R^4 + 30R^2 - 4)RCOSA
-- : [23] (35R^6 - 60R^4 + 30R^2 - 4)RSINA
2.0654e-08: [24] 70R^8 - 140R^6 + 90R^4 - 20R^2 + 1
-- : [25] R^5COS5A
-- : [26] R^5SIN5A
-- : [27] (6R^2 - 5)R^4COS4A
-- : [28] (6R^2 - 5)R^4SIN4A
-- : [29] (21R^4 - 30R^2 + 10)R^3COS3A
-- : [30] (21R^4 - 30R^2 + 10)R^3SIN3A
-- : [31] (56R^6 - 105R^4 + 60R^2 - 10)R^2COS2A
-- : [32] (56R^6 - 105R^4 + 60R^2 - 10)R^2SIN2A
-- : [33] (126R^8 - 280R^6 + 210R^4 - 60R^2 + 5)RCOSA
-- : [34] (126R^8 - 280R^6 + 210R^4 - 60R^2 + 5)RSINA
-1.2840e-09: [35] 252R^10 - 630R^8 + 560R^6 - 210R^4 + 30R^2 - 1
-- : [36] 924R^12 - 2772R^10 + 3150R^8 - 1680R^6 + 420R^4 - 42R^2 + 1
RMS OPD 0.001092 ERROR --

```

**Verification: PMMA cornea for ISO-fixture, with focus in air, with Z9000 IOL, 20.0D**

\*LENS DATA

PMMA Z9000 cornea in ISO fixture

SRF	RADIUS	THICKNESS	APERTURE RADIUS	GLASS	SPE	NOTE
OBJ	--	1.0000e+20	1.0000e+14	AIR		
1	18.568300	10.000000	8.000000	PMMA_ROD	*	
2	--	3.000000	8.000000	AIR		

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

3	--	6.000000	16.000000	BK7	C
4	--	6.250000	16.000000	WATER	C
AST	--	--	2.640233	AS	WATER
6	12.154000	1.125000	2.640233	S	HRI *
7	-12.154000	8.875000	2.516362	S	WATER
8	--	6.000000	16.000000	BK7	C
9	--	2.914830	16.000000	AIR	

IMS -- -- 2.7366e-05 S

\*CONIC AND POLYNOMIAL ASPHERIC DATA

SRF	CC	AD	AE	AF	AG
1	-0.177999	5.1185e-06	3.3540e-09	--	--
6	-1.018550	-0.000899	-1.1400e-05	--	--

\*ZERNIKE ANALYSIS

WAVELENGTH 1

Positive angle A is a rotation from the +y axis toward the +x axis.

```

-5.7675e-05: [0] 1
-- : [1] RCOSA
-- : [2] RSINA
-9.0472e-05: [3] 2R^2 - 1
-- : [4] R^2COS2A
-- : [5] R^2SIN2A
-- : [6] (3R^2 - 2)RCOSA
-- : [7] (3R^2 - 2)RSINA
-3.5900e-05: [8] 6R^4 - 6R^2 + 1
-- : [9] R^3COS3A
-- : [10] R^3SIN3A
-- : [11] (4R^2 - 3)R^2COS2A
-- : [12] (4R^2 - 3)R^2SIN2A
-- : [13] (10R^4 - 12R^2 + 3)RCOSA
-- : [14] (10R^4 - 12R^2 + 3)RSINA
-3.4880e-06: [15] 20R^6 - 30R^4 + 12R^2 - 1
-- : [16] R^4COS4A
-- : [17] R^4SIN4A
-- : [18] (5R^2 - 4)R^3COS3A
-- : [19] (5R^2 - 4)R^3SIN3A
-- : [20] (15R^4 - 20R^2 + 6)R^2COS2A
-- : [21] (15R^4 - 20R^2 + 6)R^2SIN2A
-- : [22] (35R^6 - 60R^4 + 30R^2 - 4)RCOSA
-- : [23] (35R^6 - 60R^4 + 30R^2 - 4)RSINA
-4.1784e-07: [24] 70R^8 - 140R^6 + 90R^4 - 20R^2 + 1
-- : [25] R^5COS5A
-- : [26] R^5SIN5A
-- : [27] (6R^2 - 5)R^4COS4A
-- : [28] (6R^2 - 5)R^4SIN4A
-- : [29] (21R^4 - 30R^2 + 10)R^3COS3A
-- : [30] (21R^4 - 30R^2 + 10)R^3SIN3A
-- : [31] (56R^6 - 105R^4 + 60R^2 - 10)R^2COS2A
-- : [32] (56R^6 - 105R^4 + 60R^2 - 10)R^2SIN2A
-- : [33] (126R^8 - 280R^6 + 210R^4 - 60R^2 + 5)RCOSA
-- : [34] (126R^8 - 280R^6 + 210R^4 - 60R^2 + 5)RSINA
-3.3549e-08: [35] 252R^10 - 630R^8 + 560R^6 - 210R^4 + 30R^2 - 1
-1.3002e-09: [36] 924R^12 - 2772R^10 + 3150R^8 - 1680R^6 + 420R^4 - 42R^2 + 1
RMS OPD 5.4662e-05 ERROR 1.8430e-11

```

**PMMA/Z9000 cornea for PPG-fixture.**

Apdiv 97.10

Zernike in mm

Spot diagram in image space

Pupil 2.640233

\*LENS DATA

PMMA Z9000 CORNEA PPG\_fixture n=100

SRF	RADIUS	THICKNESS	APERTURE RADIUS	GLASS	SPE	NOTE
OBJ	--	1.0000e+20	1.0000e+14	AIR		
1	11.047000	4.000000	3.000000	PMMA	*	
2	1.0000e+26	0.100000	2.650289	WATER		
AST	--	--	2.640233	AS	PUPIL	

# REPORT #: Rp2691r

## Bifocal foldable lens design based on corneal wavefront aberration

4	--	1.700000	2.640233 S	WATER
5	--	24.556458 S	2.469290 S	WATER
IMS	--	-0.140000	0.014100 S	

### \*CONIC AND POLYNOMIAL ASPHERIC DATA

SRF	CC	AD	AE	AF	AG
1	0.512084	-9.0182e-06	-4.1412e-08	--	--

### \*ZERNIKE ANALYSIS

#### WAVELENGTH 1

Positive angle A is a rotation from the +y axis toward the +x axis.

```

0.000685: [0] 1
-- : [1] RCOSA
-- : [2] RSINA
0.001280: [3] 2R^2 - 1
-- : [4] R^2COS2A
-- : [5] R^2SIN2A
-- : [6] (3R^2 - 2)RCOSA
-- : [7] (3R^2 - 2)RSINA
0.000603: [8] 6R^4 - 6R^2 + 1
-- : [9] R^4COS3A
-- : [10] R^4SIN3A
-- : [11] (4R^2 - 3)R^2COS2A
-- : [12] (4R^2 - 3)R^2SIN2A
-- : [13] (10R^4 - 12R^2 + 3)RCOSA
-- : [14] (10R^4 - 12R^2 + 3)RSINA
9.0025e-06: [15] 20R^6 - 30R^4 + 12R^2 - 1
-- : [16] R^4COS4A
-- : [17] R^4SIN4A
-- : [18] (5R^2 - 4)R^3COS3A
-- : [19] (5R^2 - 4)R^3SIN3A
-- : [20] (15R^4 - 20R^2 + 6)R^2COS2A
-- : [21] (15R^4 - 20R^2 + 6)R^2SIN2A
-- : [22] (35R^6 - 60R^4 + 30R^2 - 4)RCOSA
-- : [23] (35R^6 - 60R^4 + 30R^2 - 4)RSINA
1.4315e-07: [24] 70R^8 - 140R^6 + 90R^4 - 20R^2 + 1
-- : [25] R^5COS5A
-- : [26] R^5SIN5A
-- : [27] (6R^2 - 5)R^4COS4A
-- : [28] (6R^2 - 5)R^4SIN4A
-- : [29] (21R^4 - 30R^2 + 10)R^3COS3A
-- : [30] (21R^4 - 30R^2 + 10)R^3SIN3A
-- : [31] (56R^6 - 105R^4 + 60R^2 - 10)R^2COS2A
-- : [32] (56R^6 - 105R^4 + 60R^2 - 10)R^2SIN2A
-- : [33] (126R^8 - 280R^6 + 210R^4 - 60R^2 + 5)RCOSA
-- : [34] (126R^8 - 280R^6 + 210R^4 - 60R^2 + 5)RSINA
1.9291e-09: [35] 252R^10 - 630R^8 + 560R^6 - 210R^4 + 30R^2 - 1
-- : [36] 924R^12 - 2772R^10 + 3150R^8 - 1680R^6 + 420R^4 - 42R^2 + 1
RMS OPD 0.000787 ERROR 8.8271e-11

```

### PMMA/Z9000 cornea for PPG-fixture with 20.0D Z9000 lens.

Apdiv 97.10

Zernike in mm

Spot diagram in image space

IOL pupil 2.640233

### \*LENS DATA

PMMA Z9000 CORNEA PPG\_fixture n=100

SRF	RADIUS	THICKNESS	APERTURE RADIUS	GLASS	SPE	NOTE
OBJ	--	1.0000e+20	1.0000e+14	AIR		
1	11.047000	4.000000	3.000000	PMMA	*	
2	1.0000e+26	0.100000	2.650289 S	WATER		
AST	--	--	2.640233 AS	PUPIL		
4	--	1.700000	2.640233 S	WATER		
5	12.154000	1.103000	2.469290 S	HRI	*	
6	-12.154000	17.052946 S	3.000000	WATER		
IMS	--	-0.140000	0.019298 S			

### \*CONIC AND POLYNOMIAL ASPHERIC DATA

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

SRF	CC	AD	AZ	AP	AG
1	0.512084	-9.0182e-06	-4.1412e-08	--	--
5	-1.018548	-0.000899	-1.1412e-05	--	--

\*ZERNIKE ANALYSIS

WAVELENGTH 1

Positive angle A is a rotation from the +y axis toward the +x axis.

```

-0.000720: [0] 1
-- : [1] RCOSA
-- : [2] RSINA
-0.000632: [3] 2R^2 - 1
-- : [4] R^2COS2A
-- : [5] R^2SIN2A
-- : [6] (3R^2 - 2)RCOSA
-- : [7] (3R^2 - 2)RSINA
9.2891e-05: [8] 6R^4 - 6R^2 + 1
-- : [9] R^3COS3A
-- : [10] R^3SIN3A
-- : [11] (4R^2 - 3)R^2COS2A
-- : [12] (4R^2 - 3)R^2SIN2A
-- : [13] (10R^4 - 12R^2 + 3)RCOSA
-- : [14] (10R^4 - 12R^2 + 3)RSINA
4.8164e-06: [15] 20R^6 - 30R^4 + 12R^2 - 1
-- : [16] R^4COS4A
-- : [17] R^4SIN4A
-- : [18] (5R^2 - 4)R^3COS3A
-- : [19] (5R^2 - 4)R^3SIN3A
-- : [20] (15R^4 - 20R^2 + 6)R^2COS2A
-- : [21] (15R^4 - 20R^2 + 6)R^2SIN2A
-- : [22] (35R^6 - 60R^4 + 30R^2 - 4)RCOSA
-- : [23] (35R^6 - 60R^4 + 30R^2 - 4)RSINA
-9.1973e-09: [24] 70R^8 - 140R^6 + 90R^4 - 20R^2 + 1
-- : [25] R^5COS5A
-- : [26] R^5SIN5A
-- : [27] (6R^2 - 5)R^4COS4A
-- : [28] (6R^2 - 5)R^4SIN4A
-- : [29] (21R^4 - 30R^2 + 10)R^3COS3A
-- : [30] (21R^4 - 30R^2 + 10)R^3SIN3A
-- : [31] (56R^6 - 105R^4 + 60R^2 - 10)R^2COS2A
-- : [32] (56R^6 - 105R^4 + 60R^2 - 10)R^2SIN2A
-- : [33] (126R^8 - 280R^6 + 210R^4 - 60R^2 + 5)RCOSA
-- : [34] (126R^8 - 280R^6 + 210R^4 - 60R^2 + 5)RSINA
-1.1709e-08: [35] 252R^10 - 630R^8 + 560R^6 - 210R^4 + 30R^2 - 1
-4.6305e-10: [36] 924R^12 - 2772R^10 + 3150R^8 - 1680R^6 + 420R^4 - 42R^2 + 1
RMS OPD 0.000367 ERROR --

```

#### Appendix 4: Measurement procedure of the lens power in the waterbath

The lenses were measured for optical power in a water bath. The actual back focal length (BFL) was used as the measure of lens power in water. The lens is focused at maximum MTF at 50c/mm. The measurement is time consuming, due to the fluctuations of the MTF values. A manual through focus measurement was performed, so that the interpolated best focus could be determined off-line. For each value, an extensive settling time of 5 minutes or more was exercised.

Determination of the lens power in vivo consists of the following steps:

1. Through focus measurement in a waterbath
2. Determination of the actual BFL by interpolation
3. From the actual BFL, determine the effective focal length (EFL) by adding the longitudinal spherical aberration and the position of the 2<sup>nd</sup> principle plane.
4. Determine the lens power in water, from the equation  $1.3406/EFL$
5. Determine the lens power in vivo, from the theoretical relationship between lens power in water and lens power in vivo.

##### 1. Through focus measurement in a waterbath

A manual through focus MTF at 50c/mm was measured in a waterbath. Characteristics of the measurement set-up:

- Aperture 3mm
- Water with PEG200, refractive index 1.3406
- Slit width 10 micrometer
- 10X objective lens
- 1 through focus measurement per lens
- 5 to 8 focus positions for each focal point

##### 2. Determination of the actual BFL by interpolation

In order to be able to make an interpolation, a curve was fit through the through-focus response. For this purpose, a square double log function was used:

$$\log(MTF) = b_0 + b_1 * (\log(BFL)) + b_2 * (\log(BFL))^2$$

This function is bell-shaped by nature, similar to a normal through-focus response. The maximum of the curve is determined by  $\log(BFL) = -1/2(b_1/b_2)$

##### 3. From the actual BFL, determine the effective focal length.

In order to find the EFL, the spherical aberration has to be added to the BFL, as well as the distance between the 2<sup>nd</sup> principle plane and the posterior lens surface.

The longitudinal spherical aberration (LSA) of a bifocal lens in a waterbath is assumed identical to its monofocal equivalent. This means that the LSA can be calculated with OSLO, using the same lens geometry without the diffraction pattern, and using the *appropriate refractive indices as applicable under the measuring conditions. It means that the monofocal Z9000 lens design can be used in the calculations (only this prototype lens was measured in the waterbath).* Figure 4.1. shows the OSLO results for a 20.0D lens and the results for the 3 measured lens powers are summarised in figure 4-2.

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

*LENS DATA						
Design4ZMO 200						
SRF	RADIUS	THICKNESS	APERTURE RADIUS	GLASS	SPE	NOTE
OBJ	--	1.0000e+20	1.0000e+14	WATERMO		
AST	12.154000	1.125070	1.500000 AS	HRISILMO	*	
2	-12.154000	65.748017 S	3.000000	WATERMO		
IMS	--	--	V 6.6265e-05 S			
*CONIC AND POLYNOMIAL ASPHERIC DATA						
SRF	CC	AD	AE	AF	AG	
1	-1.018550	-0.000899	-1.1418e-05	--	--	
*WAVELENGTHS						
CURRENT	WV1/WW1	WV2/WW2	WV3/WW3			
1	0.550000	0.550000	0.550000			
	1.000000	1.000000	1.000000			
*REFRACTIVE INDICES						
SRF	GLASS	RN1	RN2	RN3	VNBR	TCE
0	WATERMO	1.340563	1.340563	1.340563	--	--
1	HRISILMO	1.463984	1.463984	1.463984	--	236.000000
2	WATERMO	1.340563	1.340563	1.340563	--	236.000000
3	IMAGE SURFACE					
minimum focus position: --						
maximum focus position: 5.000000						
number of steps: 500						
a = max MTF(f=50): 0.3557						
b = at defocus: 1.4000 mm						

Figure 4-1. OSLO output for a 20.0D lens of model Z9000, in the waterbath.

Nominal Power	Anterior radius	Posterior radius	Central Thickness	Total Power	Anterior Power	Spherical aberration
15	16.225	-16.225	1.027	15.17	7.61	2.75
19	12.797	-12.797	1.106	19.22	9.64	1.54
20	12.154	-12.154	1.125	20.23	10.16	1.40
24	10.117	-10.117	1.203	24.28	12.20	0.99
26	9.333	-9.333	1.242	26.30	13.22	0.86
30	8.075	-8.075	1.319	30.36	15.28	0.67
Units	mm	mm	mm	D	D	mm

Figure 4.2. Spherical aberration for different lens powers. Note that the spherical aberration is positive for all lens powers. This is opposite to the normal spherical lenses.

The position of the 2<sup>nd</sup> principle plane is calculated using the paraxial equations:

$$H'' = \frac{\eta_{med}}{\eta_{IOL}} * \frac{P_1}{P} * CT$$

H'' = distance between 2<sup>nd</sup> principle plane and the posterior lens surface

$\eta_{med}$  = refractive index of the surrounding medium

$\eta_{IOL}$  = refractive index of the lens

P<sub>1</sub> = power of the front surface of the lens

P = total power of the lens

CT = central thickness of the lens

Application for the monofocal Z9000 lens model results in values of 0.471, 0.517 and 0.572mm for 15.0D, 20.0D and 26D lenses respectively. According the definition of bifocal lens power, H'' is equal for far and near vision<sup>8</sup>.

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

4. Determine the lens power in water

The lens power in water is defined by  $\eta_{\text{med}}/\text{EFL}$ , with  $\eta_{\text{med}}=1.3406$  and EFL expressed in meter units. The result is not exactly identical as the lens power in vivo.

5. Determine the lens power in vivo

The lens power in vivo is defined by the lens geometry and the nominal refractive indices  $\eta_{\text{med}}=1.336$  and  $\eta_{\text{IOL}}=1.458$ . The lens power under MO-lab conditions is defined by the lens geometry and the nominal refractive indices  $\eta_{\text{med}}=1.3406$  and  $\eta_{\text{IOL}}=1.46398$ . Using the thick-lens-equation for both situations, the relationship between the two can be determined (figure 4-3).

				$\eta_{\text{med}}$ $\eta_{\text{IOL}}$	1.46398 1.34066	1.468 1.336
Power	R a	R p	Central Thickness		MO-lab Power T.L.E.	In-vivo Power T.L.E.
8.0	48.766	-48.766	0.83		8.08	8.00
8.0	40.632	-40.632	0.85		8.07	8.00
7.0	34.821	-34.821	0.87		7.08	7.00
6.0	30.463	-30.463	0.89		6.09	6.00
5.0	27.074	-27.074	0.91		5.10	5.00
10.0	24.362	-24.362	0.93		10.12	10.00
10.5	23.200	-23.200	0.94		10.82	10.50
11.0	22.143	-22.143	0.95		11.13	11.00
11.5	21.178	-21.178	0.96		11.83	11.50
12.0	20.264	-20.264	0.97		12.14	12.00
12.5	19.480	-19.480	0.98		12.84	12.50
13.0	18.729	-18.729	0.99		13.15	13.00
13.5	18.033	-18.033	1.00		13.66	13.50
14.0	17.387	-17.387	1.01		14.16	14.00
14.5	16.786	-16.786	1.02		14.67	14.50
15.0	16.225	-16.225	1.03		15.17	15.00
15.5	15.699	-15.699	1.04		15.68	15.50
16.0	15.207	-15.207	1.05		16.19	16.00
16.5	14.745	-14.745	1.06		16.69	16.50
17.0	14.309	-14.309	1.07		17.20	17.00
17.5	13.899	-13.899	1.08		17.70	17.50
18.0	13.511	-13.511	1.09		18.21	18.00
18.5	13.144	-13.144	1.10		18.71	18.50
19.0	12.797	-12.797	1.11		19.22	19.00
19.5	12.467	-12.467	1.12		19.73	19.50
20.0	12.154	-12.154	1.13		20.23	20.00
20.5	11.856	-11.856	1.13		20.74	20.50
21.0	11.572	-11.572	1.14		21.24	21.00
21.5	11.301	-11.301	1.15		21.75	21.50
22.0	11.043	-11.043	1.16		22.26	22.00
22.5	10.798	-10.798	1.17		22.78	22.50
23.0	10.566	-10.566	1.18		23.27	23.00
23.5	10.334	-10.334	1.19		23.77	23.50
24.0	10.117	-10.117	1.20		24.28	24.00
24.5	9.909	-9.909	1.21		24.78	24.50
25.0	9.709	-9.709	1.22		25.29	25.00
25.5	9.517	-9.517	1.23		25.79	25.50
26.0	9.333	-9.333	1.24		26.30	26.00
26.5	9.155	-9.155	1.25		26.81	26.50
27.0	8.984	-8.984	1.26		27.31	27.00
27.5	8.819	-8.819	1.27		27.82	27.50
28.0	8.659	-8.659	1.28		28.33	28.00
28.5	8.505	-8.505	1.29		28.84	28.51
29.0	8.357	-8.357	1.30		29.34	29.01
29.5	8.214	-8.214	1.31		29.85	29.51
30.0	8.076	-8.076	1.32		30.36	30.01

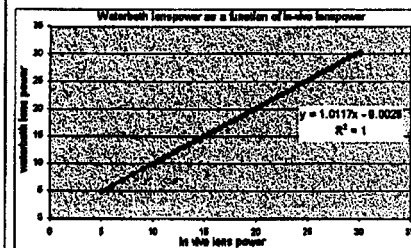


Figure 4-3. Relationship between the lens power under MO-lab conditions (waterbath) and in vivo lens power. T.L.E = thick-lens-equation.



REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

**Appendix 5: Test of the light distribution with different step heights.**

In order to determine the optimal stepheight of the diffractive profile, 3 different stepheights were tested of a 20.0D bifocal diffractive design. The target is a 50:50 light distribution.

Lens power	20.0 D	20.0 D	20.0 D
Anterior radius	12.154	12.154	12.154
Posterior radius	-12.586	-12.586	-12.586
Central thickness	1.125	1.125	1.125
R <sub>1</sub>	0.513	0.513	0.513
k	0.7065	0.6754	0.6443
cc	-1.01855	-1.01855	-1.01855
ad	0.0010482	0.0010482	0.0010482
ae	-5.7596E-06	-5.7596E-06	-5.7596E-06
stepheight (µm)	2.50	2.39	2.28

Figure 5-1. Prototype lens designs with different stepheights. These lenses have an aspherical posterior surface to correct for spherical aberration.

The lenses were tested in the Z9000 eye model, using an aperture of 3mm and water with a refractive index of 1.3405. 5 Lenses per type were tested (MO-lab service request F030). The image was focused at a maximum MTF at 50c/mm. The 4 central pixels with maximum intensity were used as a measure of the amount of light in each focus. The ratio of near and far focus determines the light distribution.

k	0.7065	0.6754	0.6443
stepheight	2.50	2.39	2.28
Lens 1	44.7%	49.9%	50.9%
Lens 2	44.8%	48.6%	50.2%
Lens 3	45.1%	48.9%	51.1%
Lens 4	46.1%	47.9%	50.2%
Lens 5	47.8%	48.2%	50.2%
Average	45.7%	48.7%	50.5%
Stdev	1.3%	0.8%	0.5%

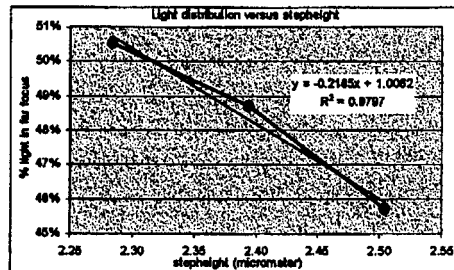


Figure 5-2. Measured percent of light in the FAR-focus versus stepheight.

The measurements show a good linearity between stepheight and light distribution. According to this data, the optimal stepheight for a 50:50% light distribution is 2.32 micrometer, which corresponds to a k-value of 0.6536.

The MTF distribution showed to be closely related to the light distribution, though the two are not identical. See figure 5-3.

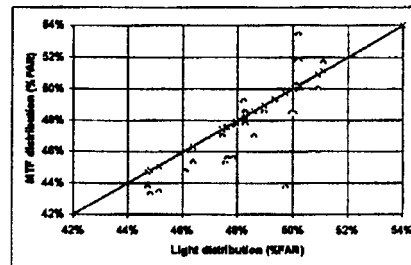


Figure 5-3. MTF(50c/mm) distribution versus light distribution.

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

### Appendix 6: Refractive indices

The performance of diffractive lenses is highly dependent on the refractive indices of the lens material and the surrounding medium and especially the difference between these two. Therefore, special care is taken to use the most accurate values for the HRI material. Furthermore, the refractive index of the water, used for measurements in the MO-lab is adjusted to the most appropriate value.

For the calculation of the image quality in vivo, the generally accepted refractive index of aqueous of 1.336 is applied. For the silicone HRI material, the measured values of the AR lab were taken as the reference ('Dispersion of HRI material', H. Weeber, PG report 2390). This is based on measurements of several batches of HRI material and it includes the dispersion of the material. The result is a refractive index of 1.45942 at 550nm and 35°C.

For the MO-lab conditions, the refractive index at 22°C and 550nm was determined using measurement data of HRI at different temperatures. This data is available in the QA-lab for a wavelength of 589nm (see figure 6-1). It is assumed that this difference of refractive index is also valid for a wavelength of 550nm.

Batch	11103	11102	11101	11003	11002	11001	11006	11005	11004	10803	10802	10801	10706	10705	10704	10703	10702	10701	Average
date	7-2-02	7-2-02	7-2-02	30-1-02	30-1-02	30-1-02	14-1-02	14-1-02	14-1-02	12-12-01	12-12-01	8-11-01	8-11-01	8-11-01	13-9-01	13-9-01	13-9-01	13-9-01	
n25	1.4616	1.4616	1.4616	1.4616	1.4615	1.4614	1.4615	1.4615	1.4614	1.4612	1.4614	1.4615	1.4618	1.4619	1.4615	1.4621	1.462	1.4619	1.461600
n37	1.4579	1.4581	1.458	1.4578	1.458	1.4572	1.4579	1.4579	1.4579	1.4579	1.4579	1.4573	1.4585	1.4586	1.4581	1.4583	1.4582	1.4581	1.457950

Figure 6-1. Measured refractive indices for 589nm and at temperatures of 25°C and 37°C.  
Raw data available at the QA lab.

Applying this to the in vivo index, the refractive index of HRI material under MO-lab conditions is 1.46398.

The refractive index of pure water, 22°C and 550nm, is 1.3342 (Handbook of Chemistry and Physics, Ed. D.R.Lide, CRC Press, Boca Raton (FL)). Under MO-lab conditions, the difference between the lens refractive index and the medium refractive index should be the same as in vivo (0.12342). Therefore, a medium refractive index of 1.34056 has to be used. For practical purposes, a tolerance of  $\pm 0.0002$  has to be added as a specification.

	Wavelength	Temperature	n <sub>lens</sub>	n <sub>medium</sub>
in vivo	550	35	1.45942	1.336
MO lab	550	22	1.46398	1.34056

Figure 6-2. Refractive indices of the lens (HRI silicone) and the medium (aqueous and water) under in vivo and MO-lab conditions, for 550nm.

In order to obtain the correct refractive index of the medium, water is mixed with 5.11% PEG200 (w/w). The mixtures used for the measurements were always measured at the AR lab on the calibrated Abbe refractometer, using a Hg light bulb with 550±10nm filter.

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

---

**Appendix 7. Raw data**

MO-lab measurement data of MTF and BFL measurements of two lens series from pilot plant service request SA020560:

K1 = bifocal foldable with optimized anterior surface

Z9000 K1 = bifocal foldable with Z9000 anterior surface

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

MECHANICAL OPTICAL LABORATORY



Pharmacia  
& Upjohn

SERVICE AANVRAAG MECHANISCH/OPTISCH LABORATORIUM

		AANVRAAGNR.: F050
AANVRAGER:	H. Weeber	AANVRAAGDATUM:
AFDELING:	AR	AFD.CODE: 570
PROJECT/OMSCHRIJVING:	Bifocal Foldable	

TYPE LENS: Prototype Bifocal Foldable

M/BATCHNUMMER	DIOPTRIE	AANTAL
SA 020560 K1	15.0D	8
SA 020560 K1	20.0D	8
SA 020560 K1	26.0D	8
SA 020560 Z9000 k1	15.0D	8
SA 020560 Z9000 k1	20.0D	8
SA 020560 Z9000 k1	26.0D	8

GEWENSTE METINGEN

MASSA		MTF/EFL	X
OVERALL DIAMETER		BACK FOCAL LENGTH	
OPTIEK DIAMETER		DIOPTRIE	
BREEDTE/DIAMETER VAN DE LUSSEN		RESOLUTIE	
DIKTE VAN DE LUSSEN		BURSTTEST	
STAND VAN DE LUSSEN/STEP HEIGHT		TREKSTERKTE	
EDGE THICKNESS			
COMPRESSION FORCE			
AXIAL DISPLACEMENT IN COMPRESSION			
TILT			
AXIAL RIGIDITY			
AMOUNT OF TOUCH			
DECENTRATION			

AFWIJKENDE METINGEN

OPMERKINGEN	Z9000 Cornea
	brekingsindex v/h water = 1.3406
	apertuur 3 en 5 mm
	focuseren bij 25 c/mm
	waarden bij 25, 50 en 100 c/mm noteren
	PH-S17-06

**REPORT #: Rp2691r**  
**Bifocal foldable lens design based on corneal wavefront aberration**

**MECHANICAL OPTICAL LABORATORY**

Servicereq. F050		Applicant: H. Weber		Dept: AR		Code: 570		Date:												
Project: 5.13		Descr: Bifocal Foldable																		
Item: SA Nummer pilot plant SA020560 (van biz-A 1m H)																				
Series: Nr. 1		No. of samples: 24		Type of dioptr: BF // 15.0/20.0/26.0 D		Remarks: Measured by MO-Laboratory														
Descr: Lense type K1 // gemeten met 3mm aperture																				
<p align="center"><b>OPTICAL BENCH</b></p> <p align="center"><b>ISO-EYE MODEL WITH Z9000 CORNEA</b></p>																				
<p align="center"><b>15.0 D</b></p>																				
<p align="center"><b>20.0 D</b></p>																				
<p align="center"><b>26.0 D</b></p>																				
<p align="center">[cm]</p>																				
unit	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
date	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
column	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
1	0.33	0.25	0.17	0.35	0.27	0.22	0.34	0.27	0.19	0.35	0.27	0.21	0.37	0.30	0.22	0.38	0.29	0.20		
2	0.37	0.30	0.22	0.36	0.29	0.23	0.34	0.26	0.22	0.33	0.26	0.20	0.36	0.28	0.21	0.40	0.30	0.23		
3	0.36	0.29	0.19	0.36	0.29	0.23	0.33	0.25	0.17	0.34	0.25	0.20	0.37	0.30	0.21	0.38	0.31	0.22		
4	0.37	0.30	0.21	0.35	0.28	0.22	0.34	0.28	0.19	0.34	0.27	0.19	0.38	0.31	0.23	0.39	0.32	0.25		
5	0.37	0.29	0.21	0.35	0.30	0.23	0.35	0.28	0.19	0.33	0.28	0.21	0.38	0.30	0.23	0.39	0.32	0.24		
6	0.35	0.29	0.17	0.35	0.28	0.21	0.34	0.30	0.22	0.35	0.28	0.23	0.36	0.30	0.22	0.37	0.30	0.22		
7	0.35	0.29	0.19	0.35	0.27	0.19	0.34	0.30	0.20	0.35	0.28	0.21	0.37	0.31	0.23	0.38	0.31	0.24		
8	0.37	0.31	0.22	0.34	0.28	0.22	0.34	0.27	0.20	0.34	0.28	0.21	0.35	0.29	0.21	0.39	0.31	0.24		
avg	0.36	0.28	0.20	0.35	0.28	0.22	0.34	0.28	0.20	0.34	0.27	0.21	0.37	0.30	0.22	0.39	0.31	0.23		
std	0.015	0.018	0.021	0.006	0.011	0.014	0.005	0.016	0.017	0.008	0.011	0.012	0.010	0.010	0.009	0.009	0.010	0.016		
min	0.33	0.25	0.17	0.34	0.27	0.19	0.33	0.25	0.17	0.33	0.25	0.19	0.35	0.28	0.21	0.37	0.29	0.20		
max	0.37	0.31	0.22	0.36	0.30	0.23	0.35	0.30	0.22	0.35	0.28	0.23	0.38	0.31	0.23	0.40	0.32	0.25		

Sign:  Check: 

**REPORT #: Rp2691r**  
**Bifocal foldable lens design based on corneal wavefront aberration**

**MECHANICAL OPTICAL LABORATORY**

Service req. F050		Applicant: H. Waaber		Dept: AR		Code: 570		Date:	
Project: 5.13		Descr: Bifocal Foldable							
Item:		SA Nummer pilot plant SA020560 (van b.z.A. tm H)							
Series: Nr. 2		No. of samples: 24		Type of dioptr? BF/15.0/20.0/26.0		Remarks: Measured by: MO-Laboratory			
Descr: Lense type K1 Z9000 // gemeten met 3mm aperture									

OPTICAL BENCH																			
ISO-EYE MODEL WITH Z9000 CORNEA																			
15.0 D						20.0 D						26.0 D							
FAR		NEAR				FAR		NEAR				FAR		NEAR					
25	50	100	25	50	100	25	50	100	25	50	100	25	50	100	25	50	100		
[c/mm]																			
1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
0.38	0.34	0.23	0.35	0.30	0.24	0.35	0.26	0.14	0.34	0.24	0.15	0.40	0.31	0.22	0.35	0.28	0.20		
0.37	0.30	0.20	0.34	0.26	0.20	0.38	0.29	0.20	0.36	0.30	0.22	0.37	0.28	0.18	0.37	0.30	0.22		
0.37	0.30	0.18	0.35	0.26	0.18	0.38	0.32	0.23	0.37	0.30	0.25	0.39	0.32	0.23	0.37	0.30	0.22		
0.36	0.31	0.21	0.36	0.30	0.24	0.40	0.31	0.22	0.38	0.30	0.21	0.39	0.33	0.25	0.36	0.31	0.25		
0.37	0.29	0.19	0.34	0.28	0.22	0.40	0.32	0.20	0.36	0.28	0.21	0.36	0.28	0.17	0.36	0.27	0.17		
0.37	0.30	0.17	0.35	0.28	0.22	0.40	0.32	0.24	0.37	0.32	0.25	0.38	0.29	0.19	0.37	0.30	0.21		
0.38	0.31	0.22	0.35	0.29	0.22	0.39	0.32	0.21	0.36	0.28	0.21	0.38	0.31	0.21	0.35	0.28	0.19		
0.37	0.30	0.21	0.34	0.28	0.22	0.38	0.29	0.18	0.36	0.30	0.24	0.38	0.31	0.21	0.38	0.30	0.22		
avg	0.37	0.31	0.20	0.35	0.28	0.22	0.39	0.30	0.36	0.29	0.22	0.38	0.30	0.21	0.36	0.29	0.21		
std	0.008	0.015	0.020	0.007	0.016	0.020	0.017	0.022	0.032	0.009	0.024	0.012	0.018	0.027	0.011	0.014	0.024		
min	0.36	0.29	0.17	0.34	0.26	0.18	0.35	0.26	0.14	0.34	0.24	0.15	0.36	0.28	0.17	0.35	0.27	0.17	
max	0.38	0.34	0.23	0.36	0.30	0.24	0.40	0.32	0.37	0.32	0.25	0.40	0.33	0.25	0.38	0.31	0.25		

page B.1 of 1

Sign:  Check: 

**REPORT #: Rp2691r**  
**Bifocal foldable lens design based on corneal wavefront aberration**

**MECHANICAL OPTICAL LABORATORY**

Servicereq.	F050	Applicant:	H. Weber	Dept:	AR	Code:	570	Date:													
Project:	5.13	Descr:	Bifocal Foldable																		
Item:		Descr:	SA Nummer pilot plant SA020560 (van bz.A /m H)																		
Series:	Nr. 1_15D	No. of samples:	24	type of diopler?	BF/15.00	Remarks:	Measured by: MO-Laboratory														
Descr:	Lentype K1 gemeten met 5mm aperture																				
total																					
8																					
OPTICAL BENCH																					
ISO-EYE MODEL WITH Z9000 CORNEA																					
FAR																					
0 degrees      -45 degrees      +45 degrees      0 degrees      -45 degrees      +45 degrees																					
25	50	100	25	50	100	25	50	100	25	50	100										
[c/mm]																					
unit																					
date																					
column	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	
1	0.40	0.35	0.27	0.40	0.33	0.23	0.42	0.34	0.23	0.31	0.27	0.18	0.33	0.25	0.14	0.31	0.28	0.21			
2	0.40	0.32	0.21	0.42	0.30	0.17	0.43	0.32	0.23	0.31	0.23	0.12	0.29	0.21	0.10	0.33	0.28	0.20			
3	0.41	0.33	0.23	0.41	0.31	0.19	0.43	0.34	0.23	0.32	0.27	0.18	0.29	0.22	0.10	0.32	0.28	0.21			
4	0.42	0.34	0.23	0.42	0.33	0.23	0.42	0.33	0.24	0.32	0.26	0.16	0.33	0.25	0.13	0.34	0.28	0.20			
5	0.41	0.32	0.20	0.38	0.27	0.15	0.39	0.30	0.21	0.30	0.24	0.13	0.28	0.20	0.09	0.30	0.25	0.18			
6	0.40	0.32	0.22	0.40	0.30	0.18	0.43	0.34	0.25	0.30	0.25	0.15	0.31	0.23	0.12	0.32	0.27	0.20			
7	0.40	0.31	0.19	0.40	0.34	0.26	0.42	0.33	0.23	0.29	0.23	0.12	0.32	0.24	0.14	0.31	0.27	0.19			
8	0.40	0.32	0.23	0.43	0.34	0.26	0.40	0.32	0.23	0.32	0.26	0.18	0.30	0.26	0.18	0.33	0.28	0.22			
avg	0.41	0.33	0.22	0.41	0.32	0.21	0.42	0.33	0.23	0.31	0.25	0.15	0.31	0.23	0.13	0.32	0.27	0.20			
std	0.008	0.013	0.024	0.016	0.024	0.042	0.015	0.014	0.011	0.011	0.016	0.027	0.019	0.021	0.029	0.013	0.011	0.012			
min	0.40	0.31	0.19	0.38	0.27	0.15	0.39	0.30	0.21	0.29	0.23	0.12	0.28	0.20	0.09	0.30	0.25	0.18			
max	0.42	0.35	0.27	0.43	0.34	0.26	0.43	0.34	0.25	0.32	0.27	0.18	0.33	0.28	0.18	0.34	0.28	0.22			
gr.avg	0.41											0.32	0.31	0.25	0.16						
gr.std	0.014											0.018		0.016	0.023	0.040					
gr.min	0.38											0.27		0.28	0.20	0.09					
gr.max	0.43											0.35		0.34	0.28	0.22					
gr.cols	(1.47)											(2.68)		(10.13.16)	(11.14.17)	(12.15.18)					

page C\_1 of 1

Sign:  Check: 

[illegible]

**Check:**



**REPORT #: Rp2691r**  
**Bifocal foldable lens design based on corneal wavefront aberration**

**MECHANICAL OPTICAL LABORATORY**

Service req.	F050	Applicant:	H. Weber	Dept:	AR	Code:	570	Date:																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																															
Project:	5.13	Descr:	Bifocal Foldable																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
Item:		Descr:	SA Nummer pilot plant SA020560 (van biz-A t/m H)																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
Series:	Nr. 1.26D	No. of samples:	24	hyper diopter?	BF// 26.00	Remarks: Measured by: MO-Laboratory																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																	
Descr:	Lentype K1 // gemeten met 5mm aperture																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																						
<table border="1"> <tr> <th colspan="2">total</th> <th colspan="18">OPTICAL BENCH</th> </tr> <tr> <th colspan="2">8</th> <th colspan="18">ISO-EYE MODEL WITH Z9000 CORNEA</th> </tr> <tr> <th colspan="2"></th> <th colspan="9">FAR</th> <th colspan="9">NEAR</th> </tr> <tr> <th colspan="2"></th> <th colspan="9">-45 degrees</th> <th colspan="9">0 degrees</th> <th colspan="9">+45 degrees</th> </tr> <tr> <th colspan="2"></th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th><th>100</th><th>25</th><th>50</th><th>100</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th><th>100</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th><th>100</th> </tr> <tr> <th>unit</th> <th>date</th> <th colspan="18">[cm]</th> </tr> <tr> <th>column</th> <th></th> <th>1</th><th>2</th><th>3</th><th>4</th><th>5</th><th>6</th><th>7</th><th>8</th><th>9</th><th>10</th><th>11</th><th>12</th><th>13</th><th>14</th><th>15</th><th>16</th><th>17</th><th>18</th><th>19</th><th>20</th> </tr> <tr> <td>1</td> <td></td> <td>0.42</td><td>0.34</td><td>0.23</td><td>0.41</td><td>0.38</td><td>0.21</td><td>0.41</td><td>0.33</td><td>0.22</td><td>0.30</td><td>0.27</td><td>0.18</td><td>0.31</td><td>0.25</td><td>0.15</td><td>0.31</td><td>0.28</td><td>0.20</td><td></td><td></td> </tr> <tr> <td>2</td> <td></td> <td>0.38</td><td>0.29</td><td>0.18</td><td>0.37</td><td>0.27</td><td>0.16</td><td>0.38</td><td>0.29</td><td>0.20</td><td>0.30</td><td>0.23</td><td>0.13</td><td>0.30</td><td>0.20</td><td>0.09</td><td>0.29</td><td>0.28</td><td>0.19</td><td></td><td></td> </tr> <tr> <td>3</td> <td></td> <td>0.40</td><td>0.29</td><td>0.17</td><td>0.41</td><td>0.33</td><td>0.22</td><td>0.38</td><td>0.29</td><td>0.18</td><td>0.29</td><td>0.25</td><td>0.15</td><td>0.31</td><td>0.25</td><td>0.17</td><td>0.29</td><td>0.24</td><td>0.16</td><td></td><td></td> </tr> <tr> <td>4</td> <td></td> <td>0.41</td><td>0.32</td><td>0.21</td><td>0.39</td><td>0.29</td><td>0.16</td><td>0.40</td><td>0.34</td><td>0.24</td><td>0.31</td><td>0.24</td><td>0.15</td><td>0.28</td><td>0.22</td><td>0.10</td><td>0.31</td><td>0.29</td><td>0.21</td><td></td><td></td> </tr> <tr> <td>5</td> <td></td> <td>0.42</td><td>0.33</td><td>0.22</td><td>0.41</td><td>0.30</td><td>0.18</td><td>0.42</td><td>0.34</td><td>0.23</td><td>0.30</td><td>0.24</td><td>0.15</td><td>0.28</td><td>0.19</td><td>0.10</td><td>0.30</td><td>0.28</td><td>0.20</td><td></td><td></td> </tr> <tr> <td>6</td> <td></td> <td>0.43</td><td>0.35</td><td>0.24</td><td>0.41</td><td>0.34</td><td>0.24</td><td>0.40</td><td>0.31</td><td>0.19</td><td>0.30</td><td>0.26</td><td>0.16</td><td>0.31</td><td>0.25</td><td>0.16</td><td>0.30</td><td>0.26</td><td>0.18</td><td></td><td></td> </tr> <tr> <td>7</td> <td></td> <td>0.41</td><td>0.31</td><td>0.20</td><td>0.40</td><td>0.32</td><td>0.22</td><td>0.40</td><td>0.30</td><td>0.18</td><td>0.31</td><td>0.26</td><td>0.16</td><td>0.31</td><td>0.26</td><td>0.17</td><td>0.30</td><td>0.24</td><td>0.15</td><td></td><td></td> </tr> <tr> <td>8</td> <td></td> <td>0.41</td><td>0.31</td><td>0.19</td><td>0.42</td><td>0.35</td><td>0.24</td><td>0.40</td><td>0.31</td><td>0.19</td><td>0.31</td><td>0.28</td><td>0.20</td><td>0.30</td><td>0.26</td><td>0.18</td><td>0.31</td><td>0.26</td><td>0.19</td><td></td><td></td> </tr> <tr> <td>avg</td> <td></td> <td>0.41</td><td>0.32</td><td>0.21</td><td>0.40</td><td>0.32</td><td>0.20</td><td>0.40</td><td>0.31</td><td>0.21</td><td>0.30</td><td>0.25</td><td>0.16</td><td>0.30</td><td>0.24</td><td>0.14</td><td>0.30</td><td>0.26</td><td>0.19</td><td></td><td></td> </tr> <tr> <td>std</td> <td></td> <td>0.012</td><td>0.022</td><td>0.024</td><td>0.016</td><td>0.035</td><td>0.033</td><td>0.014</td><td>0.021</td><td>0.022</td><td>0.007</td><td>0.017</td><td>0.021</td><td>0.011</td><td>0.028</td><td>0.037</td><td>0.008</td><td>0.018</td><td>0.021</td><td></td><td></td> </tr> <tr> <td>min</td> <td></td> <td>0.39</td><td>0.29</td><td>0.17</td><td>0.37</td><td>0.27</td><td>0.16</td><td>0.38</td><td>0.29</td><td>0.18</td><td>0.28</td><td>0.23</td><td>0.13</td><td>0.28</td><td>0.19</td><td>0.08</td><td>0.29</td><td>0.24</td><td>0.15</td><td></td><td></td> </tr> <tr> <td>max</td> <td></td> <td>0.43</td><td>0.35</td><td>0.24</td><td>0.42</td><td>0.38</td><td>0.24</td><td>0.42</td><td>0.34</td><td>0.24</td><td>0.31</td><td>0.28</td><td>0.20</td><td>0.31</td><td>0.26</td><td>0.18</td><td>0.31</td><td>0.29</td><td>0.21</td><td></td><td></td> </tr> <tr> <td>gr.avg</td> <td></td> <td colspan="20">0.40</td> </tr> <tr> <td>gr.std</td> <td></td> <td colspan="20">0.014</td> </tr> <tr> <td>gr.min</td> <td></td> <td colspan="20">0.37</td> </tr> <tr> <td>gr.max</td> <td></td> <td colspan="20">0.43</td> </tr> <tr> <td>gr.col.</td> <td></td> <td colspan="20">(1.47)</td> </tr> <tr> <td></td> <td></td> <td colspan="20">(25.8)</td> </tr> <tr> <td></td> <td></td> <td colspan="20">(3.69)</td> </tr> <tr> <td></td> <td></td> <td colspan="20">(10.13.16)</td> </tr> <tr> <td></td> <td></td> <td colspan="20">(11.14.17)</td> </tr> <tr> <td></td> <td></td> <td colspan="20">(12.15.18)</td> </tr> </table>										total		OPTICAL BENCH																		8		ISO-EYE MODEL WITH Z9000 CORNEA																				FAR									NEAR											-45 degrees									0 degrees									+45 degrees											25	50	100	25	50	100	25	50	100	25	50	100	25	50	100	25	50	100	25	50	100	unit	date	[cm]																		column		1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	1		0.42	0.34	0.23	0.41	0.38	0.21	0.41	0.33	0.22	0.30	0.27	0.18	0.31	0.25	0.15	0.31	0.28	0.20			2		0.38	0.29	0.18	0.37	0.27	0.16	0.38	0.29	0.20	0.30	0.23	0.13	0.30	0.20	0.09	0.29	0.28	0.19			3		0.40	0.29	0.17	0.41	0.33	0.22	0.38	0.29	0.18	0.29	0.25	0.15	0.31	0.25	0.17	0.29	0.24	0.16			4		0.41	0.32	0.21	0.39	0.29	0.16	0.40	0.34	0.24	0.31	0.24	0.15	0.28	0.22	0.10	0.31	0.29	0.21			5		0.42	0.33	0.22	0.41	0.30	0.18	0.42	0.34	0.23	0.30	0.24	0.15	0.28	0.19	0.10	0.30	0.28	0.20			6		0.43	0.35	0.24	0.41	0.34	0.24	0.40	0.31	0.19	0.30	0.26	0.16	0.31	0.25	0.16	0.30	0.26	0.18			7		0.41	0.31	0.20	0.40	0.32	0.22	0.40	0.30	0.18	0.31	0.26	0.16	0.31	0.26	0.17	0.30	0.24	0.15			8		0.41	0.31	0.19	0.42	0.35	0.24	0.40	0.31	0.19	0.31	0.28	0.20	0.30	0.26	0.18	0.31	0.26	0.19			avg		0.41	0.32	0.21	0.40	0.32	0.20	0.40	0.31	0.21	0.30	0.25	0.16	0.30	0.24	0.14	0.30	0.26	0.19			std		0.012	0.022	0.024	0.016	0.035	0.033	0.014	0.021	0.022	0.007	0.017	0.021	0.011	0.028	0.037	0.008	0.018	0.021			min		0.39	0.29	0.17	0.37	0.27	0.16	0.38	0.29	0.18	0.28	0.23	0.13	0.28	0.19	0.08	0.29	0.24	0.15			max		0.43	0.35	0.24	0.42	0.38	0.24	0.42	0.34	0.24	0.31	0.28	0.20	0.31	0.26	0.18	0.31	0.29	0.21			gr.avg		0.40																				gr.std		0.014																				gr.min		0.37																				gr.max		0.43																				gr.col.		(1.47)																						(25.8)																						(3.69)																						(10.13.16)																						(11.14.17)																						(12.15.18)																			
total		OPTICAL BENCH																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
8		ISO-EYE MODEL WITH Z9000 CORNEA																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
		FAR									NEAR																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																												
		-45 degrees									0 degrees									+45 degrees																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																			
		25	50	100	25	50	100	25	50	100	25	50	100	25	50	100	25	50	100	25	50	100																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																	
unit	date	[cm]																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
column		1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																		
1		0.42	0.34	0.23	0.41	0.38	0.21	0.41	0.33	0.22	0.30	0.27	0.18	0.31	0.25	0.15	0.31	0.28	0.20																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
2		0.38	0.29	0.18	0.37	0.27	0.16	0.38	0.29	0.20	0.30	0.23	0.13	0.30	0.20	0.09	0.29	0.28	0.19																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
3		0.40	0.29	0.17	0.41	0.33	0.22	0.38	0.29	0.18	0.29	0.25	0.15	0.31	0.25	0.17	0.29	0.24	0.16																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
4		0.41	0.32	0.21	0.39	0.29	0.16	0.40	0.34	0.24	0.31	0.24	0.15	0.28	0.22	0.10	0.31	0.29	0.21																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
5		0.42	0.33	0.22	0.41	0.30	0.18	0.42	0.34	0.23	0.30	0.24	0.15	0.28	0.19	0.10	0.30	0.28	0.20																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
6		0.43	0.35	0.24	0.41	0.34	0.24	0.40	0.31	0.19	0.30	0.26	0.16	0.31	0.25	0.16	0.30	0.26	0.18																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
7		0.41	0.31	0.20	0.40	0.32	0.22	0.40	0.30	0.18	0.31	0.26	0.16	0.31	0.26	0.17	0.30	0.24	0.15																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
8		0.41	0.31	0.19	0.42	0.35	0.24	0.40	0.31	0.19	0.31	0.28	0.20	0.30	0.26	0.18	0.31	0.26	0.19																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
avg		0.41	0.32	0.21	0.40	0.32	0.20	0.40	0.31	0.21	0.30	0.25	0.16	0.30	0.24	0.14	0.30	0.26	0.19																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
std		0.012	0.022	0.024	0.016	0.035	0.033	0.014	0.021	0.022	0.007	0.017	0.021	0.011	0.028	0.037	0.008	0.018	0.021																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
min		0.39	0.29	0.17	0.37	0.27	0.16	0.38	0.29	0.18	0.28	0.23	0.13	0.28	0.19	0.08	0.29	0.24	0.15																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
max		0.43	0.35	0.24	0.42	0.38	0.24	0.42	0.34	0.24	0.31	0.28	0.20	0.31	0.26	0.18	0.31	0.29	0.21																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																				
gr.avg		0.40																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
gr.std		0.014																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
gr.min		0.37																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
gr.max		0.43																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
gr.col.		(1.47)																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
		(25.8)																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
		(3.69)																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
		(10.13.16)																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
		(11.14.17)																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					
		(12.15.18)																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																																					

page E\_1 of 1

Sign. 

Check: 

**REPORT #: Rp2691r**  
**Bifocal foldable lens design based on corneal wavefront aberration**

**MECHANICAL OPTICAL LABORATORY**

Servicereq. F050		Applicant: H. Wieser		Dept: AR		Code: 570		Date:																																																																																																																																																																																
Project: 5.13		Descr: Bifocal Foldable		SA Nummer plot plant SA020560 (van blz.A t/m H)		Remarks: Measured by: MO-Laboratory																																																																																																																																																																																		
Item:		Nr. 2_15D		No. of samples: 24		Type of dioptr? BFI/15.00																																																																																																																																																																																		
Series:		Descr: Lens type Z8000 // gemeten met 5mm aperture																																																																																																																																																																																						
<p align="center"><b>OPTICAL BENCH</b></p> <p align="center"><b>ISO-EYE MODEL WITH Z8000 CORNEA</b></p>																																																																																																																																																																																								
<p align="center">[c/mm]</p>																																																																																																																																																																																								
<table border="1"> <thead> <tr> <th colspan="5">FAR</th> <th colspan="5">NEAR</th> </tr> <tr> <th colspan="5">0 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th colspan="5">-45 degrees</th> <th colspan="5">-45 degrees</th> </tr> <tr> <th colspan="5">+45 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> </tr> </thead> <tbody> <tr> <td>1</td><td>0.39</td><td>0.30</td><td>0.18</td><td>0.39</td> <td>0.41</td><td>0.32</td><td>0.20</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.31</td><td>0.27</td><td>0.19</td><td>0.31</td> </tr> <tr> <td>2</td><td>0.39</td><td>0.31</td><td>0.18</td><td>0.41</td> <td>0.33</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.17</td><td>0.33</td><td>0.29</td><td>0.23</td><td>0.32</td> </tr> <tr> <td>3</td><td>0.40</td><td>0.33</td><td>0.21</td><td>0.40</td> <td>0.32</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.21</td><td>0.32</td><td>0.29</td><td>0.22</td><td>0.32</td> </tr> <tr> <td>4</td><td>0.40</td><td>0.32</td><td>0.19</td><td>0.42</td> <td>0.36</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.29</td> <td>0.24</td><td>0.33</td><td>0.28</td><td>0.19</td><td>0.33</td> </tr> <tr> <td>5</td><td>0.40</td><td>0.32</td><td>0.22</td><td>0.42</td> <td>0.35</td><td>0.23</td><td>0.19</td><td>0.31</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.28</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>6</td><td>0.41</td><td>0.33</td><td>0.20</td><td>0.41</td> <td>0.35</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.27</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>7</td><td>0.41</td><td>0.31</td><td>0.17</td><td>0.42</td> <td>0.33</td><td>0.23</td><td>0.18</td><td>0.31</td><td>0.27</td> <td>0.19</td><td>0.31</td><td>0.25</td><td>0.16</td><td>0.32</td> </tr> <tr> <td>8</td><td>0.41</td><td>0.33</td><td>0.18</td><td>0.43</td> <td>0.34</td><td>0.22</td><td>0.22</td><td>0.32</td><td>0.26</td> <td>0.16</td><td>0.33</td><td>0.27</td><td>0.17</td><td>0.34</td> </tr> </tbody> </table>										FAR					NEAR					0 degrees					+45 degrees					-45 degrees					-45 degrees					+45 degrees					+45 degrees					25	50	100	25	50	25	50	100	25	50	25	50	100	25	50	1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31	2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32	3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32	4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33	5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32	6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32	7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32	8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34
FAR					NEAR																																																																																																																																																																																			
0 degrees					+45 degrees																																																																																																																																																																																			
-45 degrees					-45 degrees																																																																																																																																																																																			
+45 degrees					+45 degrees																																																																																																																																																																																			
25	50	100	25	50	25	50	100	25	50	25	50	100	25	50																																																																																																																																																																										
1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31																																																																																																																																																																										
2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32																																																																																																																																																																										
3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32																																																																																																																																																																										
4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33																																																																																																																																																																										
5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32																																																																																																																																																																										
6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32																																																																																																																																																																										
7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32																																																																																																																																																																										
8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34																																																																																																																																																																										
<table border="1"> <thead> <tr> <th colspan="5">FAR</th> <th colspan="5">NEAR</th> </tr> <tr> <th colspan="5">0 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th colspan="5">-45 degrees</th> <th colspan="5">-45 degrees</th> </tr> <tr> <th colspan="5">+45 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> </tr> </thead> <tbody> <tr> <td>1</td><td>0.39</td><td>0.30</td><td>0.18</td><td>0.39</td> <td>0.41</td><td>0.32</td><td>0.20</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.31</td><td>0.27</td><td>0.19</td><td>0.31</td> </tr> <tr> <td>2</td><td>0.39</td><td>0.31</td><td>0.18</td><td>0.41</td> <td>0.33</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.17</td><td>0.33</td><td>0.29</td><td>0.23</td><td>0.32</td> </tr> <tr> <td>3</td><td>0.40</td><td>0.33</td><td>0.21</td><td>0.40</td> <td>0.32</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.21</td><td>0.32</td><td>0.29</td><td>0.22</td><td>0.32</td> </tr> <tr> <td>4</td><td>0.40</td><td>0.32</td><td>0.19</td><td>0.42</td> <td>0.36</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.29</td> <td>0.24</td><td>0.33</td><td>0.28</td><td>0.19</td><td>0.33</td> </tr> <tr> <td>5</td><td>0.40</td><td>0.32</td><td>0.22</td><td>0.42</td> <td>0.35</td><td>0.23</td><td>0.19</td><td>0.31</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.28</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>6</td><td>0.41</td><td>0.33</td><td>0.20</td><td>0.41</td> <td>0.35</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.27</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>7</td><td>0.41</td><td>0.31</td><td>0.17</td><td>0.42</td> <td>0.33</td><td>0.23</td><td>0.18</td><td>0.31</td><td>0.27</td> <td>0.19</td><td>0.31</td><td>0.25</td><td>0.16</td><td>0.32</td> </tr> <tr> <td>8</td><td>0.41</td><td>0.33</td><td>0.18</td><td>0.43</td> <td>0.34</td><td>0.22</td><td>0.22</td><td>0.32</td><td>0.26</td> <td>0.16</td><td>0.33</td><td>0.27</td><td>0.17</td><td>0.34</td> </tr> </tbody> </table>										FAR					NEAR					0 degrees					+45 degrees					-45 degrees					-45 degrees					+45 degrees					+45 degrees					25	50	100	25	50	25	50	100	25	50	25	50	100	25	50	1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31	2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32	3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32	4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33	5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32	6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32	7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32	8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34
FAR					NEAR																																																																																																																																																																																			
0 degrees					+45 degrees																																																																																																																																																																																			
-45 degrees					-45 degrees																																																																																																																																																																																			
+45 degrees					+45 degrees																																																																																																																																																																																			
25	50	100	25	50	25	50	100	25	50	25	50	100	25	50																																																																																																																																																																										
1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31																																																																																																																																																																										
2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32																																																																																																																																																																										
3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32																																																																																																																																																																										
4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33																																																																																																																																																																										
5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32																																																																																																																																																																										
6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32																																																																																																																																																																										
7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32																																																																																																																																																																										
8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34																																																																																																																																																																										
<table border="1"> <thead> <tr> <th colspan="5">FAR</th> <th colspan="5">NEAR</th> </tr> <tr> <th colspan="5">0 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th colspan="5">-45 degrees</th> <th colspan="5">-45 degrees</th> </tr> <tr> <th colspan="5">+45 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> </tr> </thead> <tbody> <tr> <td>1</td><td>0.39</td><td>0.30</td><td>0.18</td><td>0.39</td> <td>0.41</td><td>0.32</td><td>0.20</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.31</td><td>0.27</td><td>0.19</td><td>0.31</td> </tr> <tr> <td>2</td><td>0.39</td><td>0.31</td><td>0.18</td><td>0.41</td> <td>0.33</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.17</td><td>0.33</td><td>0.29</td><td>0.23</td><td>0.32</td> </tr> <tr> <td>3</td><td>0.40</td><td>0.33</td><td>0.21</td><td>0.40</td> <td>0.32</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.21</td><td>0.32</td><td>0.29</td><td>0.22</td><td>0.32</td> </tr> <tr> <td>4</td><td>0.40</td><td>0.32</td><td>0.19</td><td>0.42</td> <td>0.36</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.29</td> <td>0.24</td><td>0.33</td><td>0.28</td><td>0.19</td><td>0.33</td> </tr> <tr> <td>5</td><td>0.40</td><td>0.32</td><td>0.22</td><td>0.42</td> <td>0.35</td><td>0.23</td><td>0.19</td><td>0.31</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.28</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>6</td><td>0.41</td><td>0.33</td><td>0.20</td><td>0.41</td> <td>0.35</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.27</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>7</td><td>0.41</td><td>0.31</td><td>0.17</td><td>0.42</td> <td>0.33</td><td>0.23</td><td>0.18</td><td>0.31</td><td>0.27</td> <td>0.19</td><td>0.31</td><td>0.25</td><td>0.16</td><td>0.32</td> </tr> <tr> <td>8</td><td>0.41</td><td>0.33</td><td>0.18</td><td>0.43</td> <td>0.34</td><td>0.22</td><td>0.22</td><td>0.32</td><td>0.26</td> <td>0.16</td><td>0.33</td><td>0.27</td><td>0.17</td><td>0.34</td> </tr> </tbody> </table>										FAR					NEAR					0 degrees					+45 degrees					-45 degrees					-45 degrees					+45 degrees					+45 degrees					25	50	100	25	50	25	50	100	25	50	25	50	100	25	50	1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31	2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32	3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32	4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33	5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32	6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32	7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32	8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34
FAR					NEAR																																																																																																																																																																																			
0 degrees					+45 degrees																																																																																																																																																																																			
-45 degrees					-45 degrees																																																																																																																																																																																			
+45 degrees					+45 degrees																																																																																																																																																																																			
25	50	100	25	50	25	50	100	25	50	25	50	100	25	50																																																																																																																																																																										
1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31																																																																																																																																																																										
2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32																																																																																																																																																																										
3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32																																																																																																																																																																										
4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33																																																																																																																																																																										
5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32																																																																																																																																																																										
6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32																																																																																																																																																																										
7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32																																																																																																																																																																										
8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34																																																																																																																																																																										
<table border="1"> <thead> <tr> <th colspan="5">FAR</th> <th colspan="5">NEAR</th> </tr> <tr> <th colspan="5">0 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th colspan="5">-45 degrees</th> <th colspan="5">-45 degrees</th> </tr> <tr> <th colspan="5">+45 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> </tr> </thead> <tbody> <tr> <td>1</td><td>0.39</td><td>0.30</td><td>0.18</td><td>0.39</td> <td>0.41</td><td>0.32</td><td>0.20</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.31</td><td>0.27</td><td>0.19</td><td>0.31</td> </tr> <tr> <td>2</td><td>0.39</td><td>0.31</td><td>0.18</td><td>0.41</td> <td>0.33</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.17</td><td>0.33</td><td>0.29</td><td>0.23</td><td>0.32</td> </tr> <tr> <td>3</td><td>0.40</td><td>0.33</td><td>0.21</td><td>0.40</td> <td>0.32</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.21</td><td>0.32</td><td>0.29</td><td>0.22</td><td>0.32</td> </tr> <tr> <td>4</td><td>0.40</td><td>0.32</td><td>0.19</td><td>0.42</td> <td>0.36</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.29</td> <td>0.24</td><td>0.33</td><td>0.28</td><td>0.19</td><td>0.33</td> </tr> <tr> <td>5</td><td>0.40</td><td>0.32</td><td>0.22</td><td>0.42</td> <td>0.35</td><td>0.23</td><td>0.19</td><td>0.31</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.28</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>6</td><td>0.41</td><td>0.33</td><td>0.20</td><td>0.41</td> <td>0.35</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.27</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>7</td><td>0.41</td><td>0.31</td><td>0.17</td><td>0.42</td> <td>0.33</td><td>0.23</td><td>0.18</td><td>0.31</td><td>0.27</td> <td>0.19</td><td>0.31</td><td>0.25</td><td>0.16</td><td>0.32</td> </tr> <tr> <td>8</td><td>0.41</td><td>0.33</td><td>0.18</td><td>0.43</td> <td>0.34</td><td>0.22</td><td>0.22</td><td>0.32</td><td>0.26</td> <td>0.16</td><td>0.33</td><td>0.27</td><td>0.17</td><td>0.34</td> </tr> </tbody> </table>										FAR					NEAR					0 degrees					+45 degrees					-45 degrees					-45 degrees					+45 degrees					+45 degrees					25	50	100	25	50	25	50	100	25	50	25	50	100	25	50	1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31	2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32	3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32	4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33	5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32	6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32	7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32	8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34
FAR					NEAR																																																																																																																																																																																			
0 degrees					+45 degrees																																																																																																																																																																																			
-45 degrees					-45 degrees																																																																																																																																																																																			
+45 degrees					+45 degrees																																																																																																																																																																																			
25	50	100	25	50	25	50	100	25	50	25	50	100	25	50																																																																																																																																																																										
1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31																																																																																																																																																																										
2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32																																																																																																																																																																										
3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32																																																																																																																																																																										
4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33																																																																																																																																																																										
5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32																																																																																																																																																																										
6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32																																																																																																																																																																										
7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32																																																																																																																																																																										
8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34																																																																																																																																																																										
<table border="1"> <thead> <tr> <th colspan="5">FAR</th> <th colspan="5">NEAR</th> </tr> <tr> <th colspan="5">0 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th colspan="5">-45 degrees</th> <th colspan="5">-45 degrees</th> </tr> <tr> <th colspan="5">+45 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> </tr> </thead> <tbody> <tr> <td>1</td><td>0.39</td><td>0.30</td><td>0.18</td><td>0.39</td> <td>0.41</td><td>0.32</td><td>0.20</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.31</td><td>0.27</td><td>0.19</td><td>0.31</td> </tr> <tr> <td>2</td><td>0.39</td><td>0.31</td><td>0.18</td><td>0.41</td> <td>0.33</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.17</td><td>0.33</td><td>0.29</td><td>0.23</td><td>0.32</td> </tr> <tr> <td>3</td><td>0.40</td><td>0.33</td><td>0.21</td><td>0.40</td> <td>0.32</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.21</td><td>0.32</td><td>0.29</td><td>0.22</td><td>0.32</td> </tr> <tr> <td>4</td><td>0.40</td><td>0.32</td><td>0.19</td><td>0.42</td> <td>0.36</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.29</td> <td>0.24</td><td>0.33</td><td>0.28</td><td>0.19</td><td>0.33</td> </tr> <tr> <td>5</td><td>0.40</td><td>0.32</td><td>0.22</td><td>0.42</td> <td>0.35</td><td>0.23</td><td>0.19</td><td>0.31</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.28</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>6</td><td>0.41</td><td>0.33</td><td>0.20</td><td>0.41</td> <td>0.35</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.27</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>7</td><td>0.41</td><td>0.31</td><td>0.17</td><td>0.42</td> <td>0.33</td><td>0.23</td><td>0.18</td><td>0.31</td><td>0.27</td> <td>0.19</td><td>0.31</td><td>0.25</td><td>0.16</td><td>0.32</td> </tr> <tr> <td>8</td><td>0.41</td><td>0.33</td><td>0.18</td><td>0.43</td> <td>0.34</td><td>0.22</td><td>0.22</td><td>0.32</td><td>0.26</td> <td>0.16</td><td>0.33</td><td>0.27</td><td>0.17</td><td>0.34</td> </tr> </tbody> </table>										FAR					NEAR					0 degrees					+45 degrees					-45 degrees					-45 degrees					+45 degrees					+45 degrees					25	50	100	25	50	25	50	100	25	50	25	50	100	25	50	1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31	2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32	3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32	4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33	5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32	6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32	7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32	8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34
FAR					NEAR																																																																																																																																																																																			
0 degrees					+45 degrees																																																																																																																																																																																			
-45 degrees					-45 degrees																																																																																																																																																																																			
+45 degrees					+45 degrees																																																																																																																																																																																			
25	50	100	25	50	25	50	100	25	50	25	50	100	25	50																																																																																																																																																																										
1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31																																																																																																																																																																										
2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32																																																																																																																																																																										
3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32																																																																																																																																																																										
4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33																																																																																																																																																																										
5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32																																																																																																																																																																										
6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32																																																																																																																																																																										
7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32																																																																																																																																																																										
8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34																																																																																																																																																																										
<table border="1"> <thead> <tr> <th colspan="5">FAR</th> <th colspan="5">NEAR</th> </tr> <tr> <th colspan="5">0 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th colspan="5">-45 degrees</th> <th colspan="5">-45 degrees</th> </tr> <tr> <th colspan="5">+45 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> </tr> </thead> <tbody> <tr> <td>1</td><td>0.39</td><td>0.30</td><td>0.18</td><td>0.39</td> <td>0.41</td><td>0.32</td><td>0.20</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.31</td><td>0.27</td><td>0.19</td><td>0.31</td> </tr> <tr> <td>2</td><td>0.39</td><td>0.31</td><td>0.18</td><td>0.41</td> <td>0.33</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.17</td><td>0.33</td><td>0.29</td><td>0.23</td><td>0.32</td> </tr> <tr> <td>3</td><td>0.40</td><td>0.33</td><td>0.21</td><td>0.40</td> <td>0.32</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.21</td><td>0.32</td><td>0.29</td><td>0.22</td><td>0.32</td> </tr> <tr> <td>4</td><td>0.40</td><td>0.32</td><td>0.19</td><td>0.42</td> <td>0.36</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.29</td> <td>0.24</td><td>0.33</td><td>0.28</td><td>0.19</td><td>0.33</td> </tr> <tr> <td>5</td><td>0.40</td><td>0.32</td><td>0.22</td><td>0.42</td> <td>0.35</td><td>0.23</td><td>0.19</td><td>0.31</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.28</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>6</td><td>0.41</td><td>0.33</td><td>0.20</td><td>0.41</td> <td>0.35</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.27</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>7</td><td>0.41</td><td>0.31</td><td>0.17</td><td>0.42</td> <td>0.33</td><td>0.23</td><td>0.18</td><td>0.31</td><td>0.27</td> <td>0.19</td><td>0.31</td><td>0.25</td><td>0.16</td><td>0.32</td> </tr> <tr> <td>8</td><td>0.41</td><td>0.33</td><td>0.18</td><td>0.43</td> <td>0.34</td><td>0.22</td><td>0.22</td><td>0.32</td><td>0.26</td> <td>0.16</td><td>0.33</td><td>0.27</td><td>0.17</td><td>0.34</td> </tr> </tbody> </table>										FAR					NEAR					0 degrees					+45 degrees					-45 degrees					-45 degrees					+45 degrees					+45 degrees					25	50	100	25	50	25	50	100	25	50	25	50	100	25	50	1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31	2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32	3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32	4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33	5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32	6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32	7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32	8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34
FAR					NEAR																																																																																																																																																																																			
0 degrees					+45 degrees																																																																																																																																																																																			
-45 degrees					-45 degrees																																																																																																																																																																																			
+45 degrees					+45 degrees																																																																																																																																																																																			
25	50	100	25	50	25	50	100	25	50	25	50	100	25	50																																																																																																																																																																										
1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31																																																																																																																																																																										
2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32																																																																																																																																																																										
3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32																																																																																																																																																																										
4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33																																																																																																																																																																										
5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32																																																																																																																																																																										
6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32																																																																																																																																																																										
7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32																																																																																																																																																																										
8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34																																																																																																																																																																										
<table border="1"> <thead> <tr> <th colspan="5">FAR</th> <th colspan="5">NEAR</th> </tr> <tr> <th colspan="5">0 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th colspan="5">-45 degrees</th> <th colspan="5">-45 degrees</th> </tr> <tr> <th colspan="5">+45 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> </tr> </thead> <tbody> <tr> <td>1</td><td>0.39</td><td>0.30</td><td>0.18</td><td>0.39</td> <td>0.41</td><td>0.32</td><td>0.20</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.31</td><td>0.27</td><td>0.19</td><td>0.31</td> </tr> <tr> <td>2</td><td>0.39</td><td>0.31</td><td>0.18</td><td>0.41</td> <td>0.33</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.17</td><td>0.33</td><td>0.29</td><td>0.23</td><td>0.32</td> </tr> <tr> <td>3</td><td>0.40</td><td>0.33</td><td>0.21</td><td>0.40</td> <td>0.32</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.21</td><td>0.32</td><td>0.29</td><td>0.22</td><td>0.32</td> </tr> <tr> <td>4</td><td>0.40</td><td>0.32</td><td>0.19</td><td>0.42</td> <td>0.36</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.29</td> <td>0.24</td><td>0.33</td><td>0.28</td><td>0.19</td><td>0.33</td> </tr> <tr> <td>5</td><td>0.40</td><td>0.32</td><td>0.22</td><td>0.42</td> <td>0.35</td><td>0.23</td><td>0.19</td><td>0.31</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.28</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>6</td><td>0.41</td><td>0.33</td><td>0.20</td><td>0.41</td> <td>0.35</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.27</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>7</td><td>0.41</td><td>0.31</td><td>0.17</td><td>0.42</td> <td>0.33</td><td>0.23</td><td>0.18</td><td>0.31</td><td>0.27</td> <td>0.19</td><td>0.31</td><td>0.25</td><td>0.16</td><td>0.32</td> </tr> <tr> <td>8</td><td>0.41</td><td>0.33</td><td>0.18</td><td>0.43</td> <td>0.34</td><td>0.22</td><td>0.22</td><td>0.32</td><td>0.26</td> <td>0.16</td><td>0.33</td><td>0.27</td><td>0.17</td><td>0.34</td> </tr> </tbody> </table>										FAR					NEAR					0 degrees					+45 degrees					-45 degrees					-45 degrees					+45 degrees					+45 degrees					25	50	100	25	50	25	50	100	25	50	25	50	100	25	50	1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31	2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32	3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32	4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33	5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32	6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32	7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32	8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34
FAR					NEAR																																																																																																																																																																																			
0 degrees					+45 degrees																																																																																																																																																																																			
-45 degrees					-45 degrees																																																																																																																																																																																			
+45 degrees					+45 degrees																																																																																																																																																																																			
25	50	100	25	50	25	50	100	25	50	25	50	100	25	50																																																																																																																																																																										
1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31																																																																																																																																																																										
2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32																																																																																																																																																																										
3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32																																																																																																																																																																										
4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33																																																																																																																																																																										
5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32																																																																																																																																																																										
6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32																																																																																																																																																																										
7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32																																																																																																																																																																										
8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34																																																																																																																																																																										
<table border="1"> <thead> <tr> <th colspan="5">FAR</th> <th colspan="5">NEAR</th> </tr> <tr> <th colspan="5">0 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th colspan="5">-45 degrees</th> <th colspan="5">-45 degrees</th> </tr> <tr> <th colspan="5">+45 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> </tr> </thead> <tbody> <tr> <td>1</td><td>0.39</td><td>0.30</td><td>0.18</td><td>0.39</td> <td>0.41</td><td>0.32</td><td>0.20</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.31</td><td>0.27</td><td>0.19</td><td>0.31</td> </tr> <tr> <td>2</td><td>0.39</td><td>0.31</td><td>0.18</td><td>0.41</td> <td>0.33</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.17</td><td>0.33</td><td>0.29</td><td>0.23</td><td>0.32</td> </tr> <tr> <td>3</td><td>0.40</td><td>0.33</td><td>0.21</td><td>0.40</td> <td>0.32</td><td>0.22</td><td>0.16</td><td>0.32</td><td>0.27</td> <td>0.21</td><td>0.32</td><td>0.29</td><td>0.22</td><td>0.32</td> </tr> <tr> <td>4</td><td>0.40</td><td>0.32</td><td>0.19</td><td>0.42</td> <td>0.36</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.29</td> <td>0.24</td><td>0.33</td><td>0.28</td><td>0.19</td><td>0.33</td> </tr> <tr> <td>5</td><td>0.40</td><td>0.32</td><td>0.22</td><td>0.42</td> <td>0.35</td><td>0.23</td><td>0.19</td><td>0.31</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.28</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>6</td><td>0.41</td><td>0.33</td><td>0.20</td><td>0.41</td> <td>0.35</td><td>0.25</td><td>0.25</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.33</td><td>0.27</td><td>0.20</td><td>0.32</td> </tr> <tr> <td>7</td><td>0.41</td><td>0.31</td><td>0.17</td><td>0.42</td> <td>0.33</td><td>0.23</td><td>0.18</td><td>0.31</td><td>0.27</td> <td>0.19</td><td>0.31</td><td>0.25</td><td>0.16</td><td>0.32</td> </tr> <tr> <td>8</td><td>0.41</td><td>0.33</td><td>0.18</td><td>0.43</td> <td>0.34</td><td>0.22</td><td>0.22</td><td>0.32</td><td>0.26</td> <td>0.16</td><td>0.33</td><td>0.27</td><td>0.17</td><td>0.34</td> </tr> </tbody> </table>										FAR					NEAR					0 degrees					+45 degrees					-45 degrees					-45 degrees					+45 degrees					+45 degrees					25	50	100	25	50	25	50	100	25	50	25	50	100	25	50	1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31	2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32	3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32	4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33	5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32	6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32	7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32	8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34
FAR					NEAR																																																																																																																																																																																			
0 degrees					+45 degrees																																																																																																																																																																																			
-45 degrees					-45 degrees																																																																																																																																																																																			
+45 degrees					+45 degrees																																																																																																																																																																																			
25	50	100	25	50	25	50	100	25	50	25	50	100	25	50																																																																																																																																																																										
1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31																																																																																																																																																																										
2	0.39	0.31	0.18	0.41	0.33	0.22	0.16	0.32	0.27	0.17	0.33	0.29	0.23	0.32																																																																																																																																																																										
3	0.40	0.33	0.21	0.40	0.32	0.22	0.16	0.32	0.27	0.21	0.32	0.29	0.22	0.32																																																																																																																																																																										
4	0.40	0.32	0.19	0.42	0.36	0.25	0.25	0.32	0.29	0.24	0.33	0.28	0.19	0.33																																																																																																																																																																										
5	0.40	0.32	0.22	0.42	0.35	0.23	0.19	0.31	0.27	0.18	0.33	0.28	0.20	0.32																																																																																																																																																																										
6	0.41	0.33	0.20	0.41	0.35	0.25	0.25	0.32	0.27	0.18	0.33	0.27	0.20	0.32																																																																																																																																																																										
7	0.41	0.31	0.17	0.42	0.33	0.23	0.18	0.31	0.27	0.19	0.31	0.25	0.16	0.32																																																																																																																																																																										
8	0.41	0.33	0.18	0.43	0.34	0.22	0.22	0.32	0.26	0.16	0.33	0.27	0.17	0.34																																																																																																																																																																										
<table border="1"> <thead> <tr> <th colspan="5">FAR</th> <th colspan="5">NEAR</th> </tr> <tr> <th colspan="5">0 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th colspan="5">-45 degrees</th> <th colspan="5">-45 degrees</th> </tr> <tr> <th colspan="5">+45 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th> </tr> </thead> <tbody> <tr> <td>1</td><td>0.39</td><td>0.30</td><td>0.18</td><td>0.39</td> <td>0.41</td><td>0.32</td><td>0.20</td><td>0.32</td><td>0.27</td> <td>0.18</td><td>0.31</td><td>0.27</td><td>0.19</td><td>0.31</td> </tr> <tr> <td>2</td><td>0.39</td><td>0</td></tr></tbody></table>										FAR					NEAR					0 degrees					+45 degrees					-45 degrees					-45 degrees					+45 degrees					+45 degrees					25	50	100	25	50	25	50	100	25	50	25	50	100	25	50	1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31	2	0.39	0																																																																																																						
FAR					NEAR																																																																																																																																																																																			
0 degrees					+45 degrees																																																																																																																																																																																			
-45 degrees					-45 degrees																																																																																																																																																																																			
+45 degrees					+45 degrees																																																																																																																																																																																			
25	50	100	25	50	25	50	100	25	50	25	50	100	25	50																																																																																																																																																																										
1	0.39	0.30	0.18	0.39	0.41	0.32	0.20	0.32	0.27	0.18	0.31	0.27	0.19	0.31																																																																																																																																																																										
2	0.39	0																																																																																																																																																																																						

**REPORT #: Rp2691r**  
**Bifocal foldable lens design based on corneal wavefront aberration**

**MECHANICAL OPTICAL LABORATORY**

Servicereq.	F050	Applicant:	H. Weeber	Dept:	AR	Code:	570	Date:																																																											
Project:	5.13	Descr:	Bifocal Foldable																																																																
Item:		Descr:	SA Nummer pilot plant SA020560 (van bz A.1m H)																																																																
Series:	nr_2_200	No. of samples:	24	typal diopter?	BF/20.0D	Remarks:	Measured by: MO-Laboratory																																																												
Descr:	Lentype Z9000 // gemeten met 5mm aperture																																																																		
total	8																																																																		
<b>OPTICAL BENCH</b>																																																																			
<b>ISO-EYE MODEL WITH Z9000 CORNEA</b>																																																																			
<table border="1"> <tr> <th colspan="10">FAR</th> <th colspan="10">NEAR</th> </tr> <tr> <th colspan="5">0 degrees</th> <th colspan="5">+45 degrees</th> <th colspan="5">-45 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th><th>100</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th><th>100</th> <th>25</th><th>50</th><th>100</th><th>25</th><th>50</th><th>100</th> </tr> </table>										FAR										NEAR										0 degrees					+45 degrees					-45 degrees					+45 degrees					25	50	100	25	50	100	25	50	100	25	50	100	25	50	100	25	50	100
FAR										NEAR																																																									
0 degrees					+45 degrees					-45 degrees					+45 degrees																																																				
25	50	100	25	50	100	25	50	100	25	50	100	25	50	100	25	50	100																																																		
unit	(mm)																																																																		
date																																																																			
column	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20																																															
1	0.40	0.30	0.19	0.42	0.33	0.23	0.41	0.30	0.18	0.30	0.27	0.20	0.29	0.26	0.21	0.30	0.25	0.18																																																	
2	0.40	0.31	0.18	0.42	0.32	0.20	0.40	0.28	0.15	0.30	0.27	0.18	0.31	0.27	0.17	0.31	0.25	0.15																																																	
3	0.39	0.24	0.11	0.34	0.17	0.07	0.42	0.31	0.17	0.28	0.21	0.12	0.24	0.13	0.07	0.29	0.25	0.16																																																	
4	0.40	0.30	0.16	0.41	0.31	0.17	0.40	0.28	0.14	0.31	0.26	0.17	0.30	0.25	0.15	0.30	0.25	0.16																																																	
5	0.41	0.30	0.18	0.40	0.33	0.21	0.41	0.31	0.18	0.30	0.25	0.15	0.32	0.29	0.22	0.29	0.26	0.17																																																	
6	0.38	0.28	0.17	0.41	0.28	0.14	0.42	0.32	0.21	0.30	0.27	0.18	0.29	0.24	0.14	0.32	0.28	0.23																																																	
7	0.40	0.30	0.17	0.40	0.28	0.15	0.40	0.30	0.18	0.31	0.25	0.16	0.30	0.25	0.14	0.32	0.28	0.21																																																	
8	0.40	0.31	0.19	0.41	0.32	0.21	0.40	0.30	0.18	0.31	0.28	0.19	0.30	0.26	0.16	0.30	0.26	0.18																																																	
avg	0.40	0.29	0.17	0.40	0.29	0.17	0.41	0.30	0.17	0.30	0.26	0.17	0.29	0.24	0.16	0.30	0.26	0.18																																																	
std	0.009	0.023	0.028	0.026	0.053	0.052	0.009	0.014	0.021	0.010	0.022	0.025	0.024	0.046	0.047	0.012	0.013	0.027																																																	
min	0.38	0.24	0.11	0.34	0.17	0.07	0.40	0.28	0.14	0.28	0.21	0.12	0.24	0.13	0.07	0.29	0.25	0.15																																																	
max	0.41	0.31	0.19	0.42	0.33	0.23	0.42	0.32	0.21	0.31	0.28	0.20	0.32	0.29	0.22	0.32	0.28	0.23																																																	
gr.avg	0.40			0.30			0.17			0.30			0.25			0.17																																																			
gr.std	0.016			0.033			0.034			0.016			0.031			0.034																																																			
gr.min	0.34			0.17			0.07			0.24			0.13			0.07																																																			
gr.max	0.42			0.33			0.23			0.32			0.29			0.23																																																			
gr.cols.	(1.47)			(2.58)			(3.69)			(10.13.16)			(11.14.17)			(12.15.18)																																																			

page G\_1 of 1

Sign.  Check: 

## MECHANICAL OPTICAL LABORATORY

[illegible]

page H\_-:1 of 1

**Sign.** \_\_\_\_\_

**Check:**

MECHANICAL OPTICAL LABORATORY

Service req.	F050	Applicant:	H. Weber	Dept:	AR	Code:	570	Date:	
Project:	5.13	Descr:	Bifocal Foldable						
Item:		Descr:	SA Nummer pilot plant SA020560 (van Bz. A. Um H)						
Series:	Nr. 3	No. of samples:	3	Hyper diopter?	811E/15.0/20.0/25.0	Remarks:	Measured by: MO-Laboratory		
Descr:	Lens type 811E // gemeten met 3mm apertuur								
total	OPTICAL BENCH								
3	ISO-EYE MODEL WITH Z9000 CORNEA								
	POWER		FAR		NEAR				
	25	50	100	25	50	100			
unit	(c/mm)								
date									
column	1	2	3	4	5	6	7	8	9
1	15.0 D	0.37	0.30	0.24	0.35	0.28	0.22		
2	20.0 D	0.33	0.25	0.16	0.38	0.30	0.24		
3	25.0 D	0.36	0.31	0.22	0.39	0.27	0.17		
								10	11
								12	13
								14	15
								16	17
								18	19
								20	

page 1 of 1

Sign:  Check: 

MECHANICAL OPTICAL LABORATORY

Servicereq:	F050	Applicant:	H. Weber	Dept:	AR	Code:	570	Date:																																																																													
Project:	5.13	Descr:	Bifocal Foldable																																																																																		
Item:		Descr:	SA Nimmer pilot plant SA020560 (van biz A 1m H)																																																																																		
Series:	Nr_3	No. of samples:	3	type/diopter?	811E/15.0/20.0/25.0	Remarks:	Measured by: MO-Laboratory																																																																														
Descr:	Lentype 811E // gemeten met 5mm aperture																																																																																				
Total																																																																																					
ERR																																																																																					
OPTICAL BENCH																																																																																					
ISO-EYE MODEL WITH 29000 CORNEA																																																																																					
<table border="1"> <thead> <tr> <th colspan="2"></th> <th colspan="5">FAR</th> <th colspan="5">NEAR</th> </tr> <tr> <th colspan="2"></th> <th colspan="5">-45 degrees</th> <th colspan="5">+45 degrees</th> </tr> <tr> <th colspan="2"></th> <th>0 degrees</th> <th>25</th> <th>50</th> <th>100</th> <th>100</th> <th>25</th> <th>50</th> <th>100</th> <th>25</th> <th>50</th> <th>100</th> </tr> </thead> <tbody> <tr> <td>POWER</td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> </tr> <tr> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> </tr> <tr> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> </tr> </tbody> </table>												FAR					NEAR							-45 degrees					+45 degrees							0 degrees	25	50	100	100	25	50	100	25	50	100	POWER																																						
		FAR					NEAR																																																																														
		-45 degrees					+45 degrees																																																																														
		0 degrees	25	50	100	100	25	50	100	25	50	100																																																																									
POWER																																																																																					
unit		[c/mm]																																																																																			
date																																																																																					
column		1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20																																																																
1	15.0 D	0.20	0.08	0.01		0.20	0.10	0.02	0.19	0.10	0.02	0.23	0.07	0.01	0.22	0.08	0.01	0.23	0.08	0.00																																																																	
2	20.0 D	0.20	0.10	0.02		0.19	0.09	0.02	0.19	0.09	0.02	0.22	0.07	0.04	0.23	0.09	0.02	0.22	0.08	0.02																																																																	
3	25.0 D	0.18	0.06	0.04		0.18	0.07	0.03	0.18	0.07	0.02	0.16	0.03	0.03	0.14	0.02	0.03	0.15	0.03	0.02																																																																	

Sign.  Check: 

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

MECHANICAL OPTICAL LABORATORY

Servicereq. F050	Applicant: H.Weeber	Dept: AR	Code: 570	Date: [REDACTED]
Project: 5.13	Descr: Bifocal Foldable			
Item:	Descr: SA Nummer pilot plant SA020560 (van blz.A t/m H)			

	METINGEN	Spec.Doc.	WVS
LENZEN	MTF:	geen	geen

VERBRUIKTE LENZEN

BIFOCAL FOLDABLE			
SERVICE NR	MODEL	DIOPTER	Gemeten Nummers
SA 020560 K1	513BF	15.0D	2/5 en 7/10
SA 020560 K1	513BF	20.0D	1/8
SA 020560 K1	513BF	26.0D	1/8
SA 020560 K1Z900	513BF	15.0D	1/8
SA 020560 K1Z900	513BF	20.0D	1/8
SA 020560 K1Z900	513BF	26.0D	1/8

811E		
Batch nr.	Lens nr.	Power
4224 06	01 9	15.0 D
4225 21	01 7	20.0D
4207 22	11 2	25.0 D

ALGEMEEN:

De MTF metingen zijn gedaan met:

ISO houder en Z9000 cornea  
Objectief 20.8X  
lenshouder 3-en 5mm

Gullstrand houder met Z9000 cornea  
Objectief 39.3X  
lensinserts 3- en 5mm

NUT [REDACTED] [REDACTED]

Vloeistof is speciaal aangemaakt:

millipure met PEG (Brekingsindex:  $1.3406 \pm 0.0002$ )

MEETPROGRAMM

Hetzelfde als bij MTF metingen

De window was handmatig ingesteld

Hij is bij het array signal op het midden van de laatste 2 blokjes afgesteld

DE METINGEN

Er is alleen bij 25 c/mm gefocuseerd

De getallen bij 50 en 100c/mm komen van het focuspunt 25c/mm

page Comments Sign. [REDACTED]

Check: [REDACTED]

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

MECHANICAL OPTICAL LABORATORY



Pharmacia  
& Upjohn

		AANVRAAGNR.: F054
AANVRAGER:	H. Weeber	AANVRAAGDATUM:
AFDELING:	AR	AFD.CODE:
PROJECT/OMSCHRIJVING:	5.13 Bifocal Foldable	

TYPE LENS: Prototype bifocal foldable

M-/BATCHNUMMER	DIOPTRIE	AANTAL
SA 020560 Z9000 K1	15.0D	8
SA 020560 Z9000 K1	20.0D	8
SA 020560 Z9000 K1	26.0D	8

GEWENSTE METINGEN

MASSA		MTF/EFL	
OVERALL DIAMETER		BACK FOCAL LENGTH	X
OPTIEK DIAMETER		DIOPTRIE	
BREEDTE-/DIAMETER VAN DE LUSSEN		RESOLUTIE	
DIKTE VAN DE LUSSEN		BURSTTEST	
STAND VAN DE LUSSEN/STEP HEIGHT		TREKSTERKTE	
EDGE THICKNESS			
COMPRESSION FORCE			
AXIAL DISPLACEMENT IN COMPRESSION			
TILT			
AXIAL RIGIDITY			
AMOUNT OF TOUCH			
DECENTRATION			

AFWIJKENDE METINGEN

OPMERKINGEN	Meting in de watebak
	apertuur 3 mm
	Brekkingsindex van het water = 1.3406 ±0.0002



**REPORT #: Rp2691r**  
**Bifocal foldable lens design based on corneal wavefront aberration**

**MECHANICAL OPTICAL LABORATORY**

Servicereq. F054	Applicant: H. Weeber	Dept: AR	Date: [REDACTED]
Project: 5.13	Descr: Multifocal foldable lens // Prototype K1-29000	Code: [REDACTED]	
Item: SERIE_1A	No. of samples: 4	Type: dioptr/7	Remarks: Measured by: MO laboratory
Series: K1-29000/15.00			
Descr: [REDACTED]			

EROS OPTICAL BENCH																				
WATERBATH // WITH FIXTURE APERTURE = 3mm																				
		Lens-1				Lens-2				Lens-3				Lens-4						
		NEAR		FAR		NEAR		FAR		NEAR		FAR		NEAR		FAR				
		BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF			
		50c/mm				50c/mm				50c/mm				50c/mm						
unit	date	50c/mm																		
column	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
1	75.086	0.02	95.767	0.06		75.125	0.16	96.837	0.19		76.109	0.01	98.685	0.06		76.000	0.07	96.752	0.06	
2	74.203	0.10	94.616	0.18		74.221	0.24	95.706	0.26		75.251	0.16	95.597	0.15		74.651	0.15	95.750	0.17	
3	73.400	0.29	93.505	0.17		73.333	0.41	94.635	0.27		74.302	0.35	94.488	0.20		74.200	0.22	94.750	0.20	
4	72.602	0.27	92.404	0.08		72.441	0.27	93.525	0.19		73.356	0.21	93.395	0.11		73.853	0.33	93.750	0.15	
5	71.801	0.14	91.302	0.07		71.553	0.16	92.400	0.14		72.408	0.15	92.296	0.07		73.300	0.28	92.750	0.10	
6	71.003	0.02	90.201	0.03		70.661	0.09	91.306	0.08		71.451	0.01	91.194	0.07		72.850	0.17	91.750	0.10	
7	70.198		89.104			69.774		90.215			70.502		90.168			72.200		91.750		
8	69.402		88.874			69.199		89.934			70.407					71.750		90.803		

Sign: [REDACTED] Check: [REDACTED]

page A\_1 of 1

**REPORT #: Rp2691r**  
**Bifocal foldable lens design based on corneal wavefront aberration**

**MECHANICAL OPTICAL LABORATORY**

Servicereq. F054		Applicant: H.Weeber		Dept: AR		Code:		Date:	
Project: 5.13		Descr: Multifocal foldable lens // Prototype K1-Z9000							
Item:		Descr:							
Series: SERIE_1B		No. of samples: 4		Type of dioptr: 7		K1-Z9000/15.0D		Remarks: Measured by: MO laboratory	
Descr:									
<div style="display: flex; justify-content: space-between;"> <div> <b>total</b>  <b>ERR</b> </div> <div> <b>EROS OPTICAL BENCH</b>  <b>WATERBATH // WITH FIXTURE APERTURE = 3mm</b> </div> </div>									
		Lens-5		Lens-6		Lens-7		Lens-8	
		NEAR		NEAR		NEAR		NEAR	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
unit		BFL		BFL		BFL		BFL	
date		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	
		BFL		BFL		BFL		BFL	
		MTF		MTF		MTF		MTF	
		FAR		FAR		FAR		FAR	
		50c/mm		50c/mm		50c/mm		50c/mm	

**REPORT #: Rp2691r**  
**Bifocal foldable lens design based on corneal wavefront aberration**

**MECHANICAL OPTICAL LABORATORY**

Service req:	F054	Applicant:	H. Weber	Dept:	AR	Code:		Date:	
Project:	5.13	Descr:	Multifocal foldable lens // Prototype K1-Z8000						
Item:		Descr:							
Series:	SERIE_2A	No. of samples:	4	type of dioptr?	K1-Z8000/20.00	Remarks:	Measured by: MO laboratory		
Descr:									

EROS OPTICAL BENCH																				
WATERBATH // WITH FIXTURE APERTURE = 3mm																				
		Lens-1				Lens-2				Lens-3				Lens-4						
		NEAR		FAR		NEAR		FAR		NEAR		FAR		NEAR		FAR				
		BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF			
		50c/mm				50c/mm				50c/mm				50c/mm						
unit	date																			
total	ERR																			
column	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
1	53.598	64.708	65.803	0.05		53.362	0.01	69.803	0.02		60.876	0.03	66.183	0.03		53.510	0.02	70.910	0.05	
2	54.592	0.13	66.904	0.08		54.456	0.11	68.803	0.14		59.622	0.06	67.112	0.07		54.510	0.11	66.405	0.06	
3	55.597	0.22	67.504	0.11		55.561	0.16	68.301	0.17		58.360	0.16	68.612	0.12		55.508	0.24	67.504	0.13	
4	56.097	0.24	68.052	0.14		56.106	0.21	67.805	0.18		57.154	0.29	69.110	0.18		56.008	0.32	68.107	0.20	
5	56.596	0.12	68.602	0.16		56.401	0.28	67.207	0.08		56.553	0.20	69.613	0.15		56.323	0.30	68.657	0.20	
6	57.092	0.05	69.153	0.09		56.701	0.28	64.705			55.952	0.08	70.109	0.12		56.609	0.15	69.204	0.14	
7	57.592	0.03	69.703	0.03		57.004	0.17				55.349		70.607	0.05		57.110	0.03	69.753	0.04	
8	58.596		71.248			57.310	0.11				53.326		72.137			58.106				
9	59.864					59.873										59.629				
10																				

page C.1 of 1

Sign.  Check: 

## MECHANICAL OPTICAL LABORATORY

[illegible]

page D\_:1 of 1

Sign. [Signature] Check: 99

**REPORT #: Rp2691r**  
**Bifocal foldable lens design based on corneal wavefront aberration**

**MECHANICAL OPTICAL LABORATORY**

Service req.:	F054	Applicant:	H. Weber	Dept:	AR	Code:		Date:	
Project:	5.13	Descr:	Multifocal foldable lens // Prototype K1-Z9000						
Item:		Descr:							
Series:	SERIE_3A	No. of samples:	4	Type of dioplet?	K1-Z9000/28.00	Remarks:	Measured by: MO laboratory		
Descr:									
EROS OPTICAL BENCH									
WATERBATH // WITH FIXTURE APERTURE = 3mm									
Lens-1									
NEAR		FAR		NEAR		FAR		NEAR	
BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF
500mm		500mm		500mm		500mm		500mm	
1	2	3	4	5	6	7	8	9	10
47.214	0.04	48.617	0.05		45.665	0.13	49.015	0.07	
45.804	0.08	50.615	0.12		45.415	0.18	50.784	0.13	
45.403	0.25	51.117	0.19		45.165	0.29	51.265	0.22	
45.003	0.35	51.616	0.21		44.816	0.32	51.766	0.19	
44.805	0.23	52.112	0.10		44.862	0.28	52.265	0.08	
44.205	0.09	52.617	0.03		44.415	0.19	52.765	0.11	
43.805	0.08	53.118			44.163	0.11	54.547		
43.405		54.485			43.912				
41.813					41.644				
Lens-2									
NEAR		FAR		NEAR		FAR		NEAR	
BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF
500mm		500mm		500mm		500mm		500mm	
11	12	13	14	15	16	17	18	19	20
47.179	0.01	49.469	0.04		48.036	0.06	48.286	0.07	
45.674	0.05	50.309	0.07		47.057	0.04	50.789	0.17	
45.424	0.14	50.752	0.14		46.041	0.18	51.287	0.27	
45.173	0.23	51.218	0.23		45.276	0.26	51.787	0.20	
44.924	0.31	51.653	0.21		45.022	0.37	52.286	0.12	
44.673	0.22	52.106	0.15		44.772	0.33	52.788	0.08	
44.425	0.11	52.335	0.09		44.522	0.27	53.286		
44.175		52.556			44.270	0.140	54.330		
43.922		54.934			44.023				
41.164					43.772				

page E\_1 of 1

Sign:  Check: 

**REPORT #: Rp2691r**  
**Bifocal foldable lens design based on corneal wavefront aberration**

**MECHANICAL OPTICAL LABORATORY**

Service req. F054		Applicant: H. Weber		Dept: AR		Code:		Date:																																																																																																																																																																																																																																																																																																																																																															
Project: 5.13		Descr: Multifocal foldable lens // Prototype K1-Z9000																																																																																																																																																																																																																																																																																																																																																																					
Item:		Descr:																																																																																																																																																																																																																																																																																																																																																																					
Series: SERIE3B		No. of samples: 4		Type/ diopter? 1K-Z9000/26.00		Remarks:		Measured by: MO laboratory																																																																																																																																																																																																																																																																																																																																																															
Descr:																																																																																																																																																																																																																																																																																																																																																																							
<p align="center"><b>EROS OPTICAL BENCH</b></p> <p align="center">WATERBATH // WITH FIXTURE APERTURE = 3mm</p>																																																																																																																																																																																																																																																																																																																																																																							
<table border="1"> <thead> <tr> <th colspan="2"></th> <th colspan="4">Lens-6</th> <th colspan="4">Lens-7</th> <th colspan="4">Lens-8</th> <th colspan="4">Lens-9</th> </tr> <tr> <th colspan="2"></th> <th colspan="2">NEAR</th> <th colspan="2">FAR</th> <th colspan="2">NEAR</th> <th colspan="2">FAR</th> <th colspan="2">NEAR</th> <th colspan="2">FAR</th> <th colspan="2">NEAR</th> <th colspan="2">FAR</th> </tr> <tr> <th colspan="2"></th> <th>BFL</th> <th>MTF</th> <th>BFL</th> <th>MTF</th> <th>BFL</th> <th>MTF</th> <th>BFL</th> <th>MTF</th> <th>BFL</th> <th>MTF</th> <th>BFL</th> <th>MTF</th> <th>BFL</th> <th>MTF</th> <th>BFL</th> <th>MTF</th> </tr> <tr> <th colspan="2"></th> <th colspan="4">50c/mm</th> <th colspan="4">50c/mm</th> <th colspan="4">50c/mm</th> <th colspan="4">50c/mm</th> </tr> </thead> <tbody> <tr> <td>unit</td> <td></td> <td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td> </tr> <tr> <td>date</td> <td></td> <td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td><td></td> </tr> <tr> <td>column</td> <td></td> <td>1</td><td>2</td><td>3</td><td>4</td><td>5</td><td>6</td><td>7</td><td>8</td><td>9</td><td>10</td><td>11</td><td>12</td><td>13</td><td>14</td><td>15</td><td>16</td><td>17</td><td>18</td><td>19</td><td>20</td> </tr> <tr> <td>1</td> <td></td> <td>48.081</td><td></td><td>48.696</td><td></td><td></td><td>47.988</td><td></td><td>49.287</td><td></td><td></td><td>46.981</td><td></td><td>49.361</td><td></td><td></td><td>46.922</td><td></td><td>48.252</td><td></td><td></td> </tr> <tr> <td>2</td> <td></td> <td>45.454</td><td>0.09</td><td>50.697</td><td>0.05</td><td></td><td>45.600</td><td>0.08</td><td>50.288</td><td>0.04</td><td></td><td>45.483</td><td>0.10</td><td>50.882</td><td>0.06</td><td></td><td>45.673</td><td>0.05</td><td>50.754</td><td>0.07</td><td></td> </tr> <tr> <td>3</td> <td></td> <td>45.055</td><td>0.27</td><td>51.192</td><td>0.12</td><td></td><td>45.454</td><td>0.11</td><td>50.782</td><td>0.06</td><td></td><td>45.235</td><td>0.18</td><td>51.362</td><td>0.14</td><td></td><td>45.424</td><td>0.15</td><td>51.254</td><td>0.16</td><td></td> </tr> <tr> <td>4</td> <td></td> <td>44.543</td><td>0.31</td><td>51.695</td><td>0.23</td><td></td><td>45.203</td><td>0.21</td><td>51.286</td><td>0.16</td><td></td><td>44.982</td><td>0.34</td><td>51.613</td><td>0.19</td><td></td><td>45.172</td><td>0.27</td><td>51.753</td><td>0.27</td><td></td> </tr> <tr> <td>5</td> <td></td> <td>44.025</td><td>0.14</td><td>52.195</td><td>0.23</td><td></td><td>44.952</td><td>0.34</td><td>51.784</td><td>0.27</td><td></td><td>44.732</td><td>0.39</td><td>51.862</td><td>0.27</td><td></td><td>44.984</td><td>0.37</td><td>52.003</td><td>0.28</td><td></td> </tr> <tr> <td>6</td> <td></td> <td>43.532</td><td>0.09</td><td>52.694</td><td>0.11</td><td></td><td>44.702</td><td>0.32</td><td>52.035</td><td>0.31</td><td></td><td>44.482</td><td>0.32</td><td>52.112</td><td>0.25</td><td></td><td>44.734</td><td>0.33</td><td>52.254</td><td>0.26</td><td></td> </tr> <tr> <td>7</td> <td></td> <td>43.031</td><td>0.08</td><td>54.606</td><td></td><td></td><td>44.454</td><td>0.23</td><td>52.284</td><td>0.27</td><td></td><td>44.232</td><td>0.22</td><td>52.362</td><td>0.22</td><td></td><td>44.481</td><td>0.27</td><td>52.754</td><td>0.13</td><td></td> </tr> <tr> <td>8</td> <td></td> <td>42.452</td><td></td><td></td><td></td><td></td><td>44.202</td><td>0.20</td><td>52.534</td><td>0.16</td><td></td><td>43.982</td><td>0.13</td><td>52.612</td><td>0.18</td><td></td><td>44.233</td><td>0.18</td><td>53.254</td><td>0.090</td><td></td> </tr> <tr> <td>9</td> <td></td> <td></td><td></td><td></td><td></td><td></td><td>43.953</td><td>0.12</td><td>52.784</td><td>0.09</td><td></td><td>43.730</td><td>0.07</td><td>52.863</td><td>0.06</td><td></td><td>44.983</td><td>0.12</td><td>54.755</td><td></td><td></td> </tr> <tr> <td>10</td> <td></td> <td></td><td></td><td></td><td></td><td></td><td>41.930</td><td></td><td>54.787</td><td></td><td></td><td>42.853</td><td></td><td>54.234</td><td></td><td></td><td>42.421</td><td></td><td></td><td></td><td></td> </tr> </tbody> </table>												Lens-6				Lens-7				Lens-8				Lens-9						NEAR		FAR		NEAR		FAR		NEAR		FAR		NEAR		FAR				BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF			50c/mm				50c/mm				50c/mm				50c/mm				unit																		date																		column		1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	1		48.081		48.696			47.988		49.287			46.981		49.361			46.922		48.252			2		45.454	0.09	50.697	0.05		45.600	0.08	50.288	0.04		45.483	0.10	50.882	0.06		45.673	0.05	50.754	0.07		3		45.055	0.27	51.192	0.12		45.454	0.11	50.782	0.06		45.235	0.18	51.362	0.14		45.424	0.15	51.254	0.16		4		44.543	0.31	51.695	0.23		45.203	0.21	51.286	0.16		44.982	0.34	51.613	0.19		45.172	0.27	51.753	0.27		5		44.025	0.14	52.195	0.23		44.952	0.34	51.784	0.27		44.732	0.39	51.862	0.27		44.984	0.37	52.003	0.28		6		43.532	0.09	52.694	0.11		44.702	0.32	52.035	0.31		44.482	0.32	52.112	0.25		44.734	0.33	52.254	0.26		7		43.031	0.08	54.606			44.454	0.23	52.284	0.27		44.232	0.22	52.362	0.22		44.481	0.27	52.754	0.13		8		42.452					44.202	0.20	52.534	0.16		43.982	0.13	52.612	0.18		44.233	0.18	53.254	0.090		9							43.953	0.12	52.784	0.09		43.730	0.07	52.863	0.06		44.983	0.12	54.755			10							41.930		54.787			42.853		54.234			42.421				
		Lens-6				Lens-7				Lens-8				Lens-9																																																																																																																																																																																																																																																																																																																																																									
		NEAR		FAR		NEAR		FAR		NEAR		FAR		NEAR		FAR																																																																																																																																																																																																																																																																																																																																																							
		BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF	BFL	MTF																																																																																																																																																																																																																																																																																																																																																						
		50c/mm				50c/mm				50c/mm				50c/mm																																																																																																																																																																																																																																																																																																																																																									
unit																																																																																																																																																																																																																																																																																																																																																																							
date																																																																																																																																																																																																																																																																																																																																																																							
column		1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20																																																																																																																																																																																																																																																																																																																																																		
1		48.081		48.696			47.988		49.287			46.981		49.361			46.922		48.252																																																																																																																																																																																																																																																																																																																																																				
2		45.454	0.09	50.697	0.05		45.600	0.08	50.288	0.04		45.483	0.10	50.882	0.06		45.673	0.05	50.754	0.07																																																																																																																																																																																																																																																																																																																																																			
3		45.055	0.27	51.192	0.12		45.454	0.11	50.782	0.06		45.235	0.18	51.362	0.14		45.424	0.15	51.254	0.16																																																																																																																																																																																																																																																																																																																																																			
4		44.543	0.31	51.695	0.23		45.203	0.21	51.286	0.16		44.982	0.34	51.613	0.19		45.172	0.27	51.753	0.27																																																																																																																																																																																																																																																																																																																																																			
5		44.025	0.14	52.195	0.23		44.952	0.34	51.784	0.27		44.732	0.39	51.862	0.27		44.984	0.37	52.003	0.28																																																																																																																																																																																																																																																																																																																																																			
6		43.532	0.09	52.694	0.11		44.702	0.32	52.035	0.31		44.482	0.32	52.112	0.25		44.734	0.33	52.254	0.26																																																																																																																																																																																																																																																																																																																																																			
7		43.031	0.08	54.606			44.454	0.23	52.284	0.27		44.232	0.22	52.362	0.22		44.481	0.27	52.754	0.13																																																																																																																																																																																																																																																																																																																																																			
8		42.452					44.202	0.20	52.534	0.16		43.982	0.13	52.612	0.18		44.233	0.18	53.254	0.090																																																																																																																																																																																																																																																																																																																																																			
9							43.953	0.12	52.784	0.09		43.730	0.07	52.863	0.06		44.983	0.12	54.755																																																																																																																																																																																																																																																																																																																																																				
10							41.930		54.787			42.853		54.234			42.421																																																																																																																																																																																																																																																																																																																																																						

page F.1 of 1

Check:

Sign:

REPORT #: Rp2691r  
Bifocal foldable lens design based on corneal wavefront aberration

MECHANICAL OPTICAL LABORATORY

Servicereq.	F054	Applicant:	H. Weeber	Dept:	AR	Code:		Date:	30-08-2002
Project:	5.13	Descr:	Multifocal foldable lens // Prototype K1-Z9000						
Item:		Descr:							

Verbruikte Lenzen:

Position	SA numb	Power	Lens nr.
1	SA02056	15.0D	1
2	SA02056	15.0D	2
3	SA02056	15.0D	3
4	SA02056	15.0D	4
5	SA02056	15.0D	5
6	SA02056	15.0D	6
7	SA02056	15.0D	7
8	SA02056	15.0D	8

SA numb	Power	Lens nr.
SA02056	20.0D	1
SA02056	20.0D	2
SA02056	20.0D	3
SA02056	20.0D	4
SA02056	20.0D	5
SA02056	20.0D	6
SA02056	20.0D	7
SA02056	20.0D	8

SA numb	Power	Lens nr.
SA02056	26.0D	1
SA02056	26.0D	2
SA02056	26.0D	3
SA02056	26.0D	5
SA02056	26.0D	6
SA02056	26.0D	7
SA02056	26.0D	8
SA02056	26.0D	9

Algemeen:

De Back Focal Length is als volgt gemeten.

BENODIGDHEDEN: Waterbak voor BFL  
10X objectief  
Lenshouder (dezelfde als voor de ISO metingen 3mm)

MEETPROGRAMM Hetzelfde als bij MTF metingen  
De window was handmatig ingesteld  
Hij is bij het array signal op het midden van de laatste 2 blokjes afgesteld

WATERBAK Brekingsindex:  $1.3406 \pm 0.0002$   
De vloeistof is apart voor deze lenzen aangemaakt.  
Millipure + PEG

Verdere afstellingen van de opstelling als normaal

BFL meting:  
De meting is niet op één plaats uitgevoerd maar op meerdere plaatsen  
dit is als volgt gedaan:  
Op het opp. van de lens is de lengteunit genuld.  
Daarna de waterbak zover verplaatst dat er net een MTF waarde zichtbaar werd  
Deze verplaatsing is genoteerd  
Daarna in verschillende stappen door het best focus heen totdat er geen MTF meer te zien was.  
Genoteerd wordt: de stapgrootte en de MTF waarde  
Door H. Weeber is hier een best BFL uit berekend.

page Comments

Sign.

Check:

**Appendix 8. CD rom: data files**  
(on file in at Documentation Control)